

Expander 5.0.2 Online Documentation

Table of Contents

Introduction	1
Starting EXPANDER	2
Input Data	3
Tabular Data File	3
CEL Files	5
Gene Sets	8
Preprocessing GE Data	9
Viewing Data Plots	14
Clustering GE Data	16
Biclustering GE Data	19
Network Based Grouping of GE Data	24
Group Analysis Tools	26
Functional Analysis	26
Promoter Analysis	30
Location Enrichment Analysis	34
miRNA Targets Enrichment Analysis	38
General Enrichment Analysis	42
Network Based Group Analysis	45
Matrix Visualizations (Heat maps)	46
Hierarchical Clustering and Visualization	47
PCA Transformation	49
Additional Options	49
File Formats	52
Sample Input Files	54
Supplied Files	55
Settings	61
R External Application	61
FAQ	62
Copyrights Information	66
References	67

Introduction

EXPANDER (EXpression Analyzer and DisplayER) is a java-based tool for analysis of gene expression data. It is capable of (1) preprocessing (2) visualizing (3) clustering (4) biclustering and (5) performing downstream analysis of clusters and biclusters such as functional enrichment and promoter analysis (i.e. analysis of gene groups for enrichment of transcription factor binding sites in their promoters).

EXPANDER incorporates several conventional gene expression analysis algorithms and custom ones that have been developed in the computational genomics group in Tel-Aviv University, and provides them with an easy-to-operate user interface.

EXPANDER versions are available for Windows OS and for Linux/Unix OS and require the pre-installation of the Java Runtime Environment (JRE) 5.0 (or higher). The Java Runtime Environment can be installed via: <http://java.sun.com/javase/downloads/index.jsp>.


The CEL file preprocessing and the newly added SAM filter utilities require the pre-installation of one of the recent versions of R, a free software environment for statistical computing and graphics. For installation instructions, please refer to [R External Application](#) section.

Starting EXPANDER

Double click on the **Expander.bat** file, which is located under the Expander directory (alternatively, in Linux, open a Terminal window, cd into the Expander directory, and run the command: './Expander.bat').

When running on Linux/Unix OS, make sure that you have rwx permissions for the Expander directory and for the directory in which your data is located. Also make sure that you have rx permissions for all *.exe files that are under your Expander directory.

Upon running the program, the main menu bar appears:

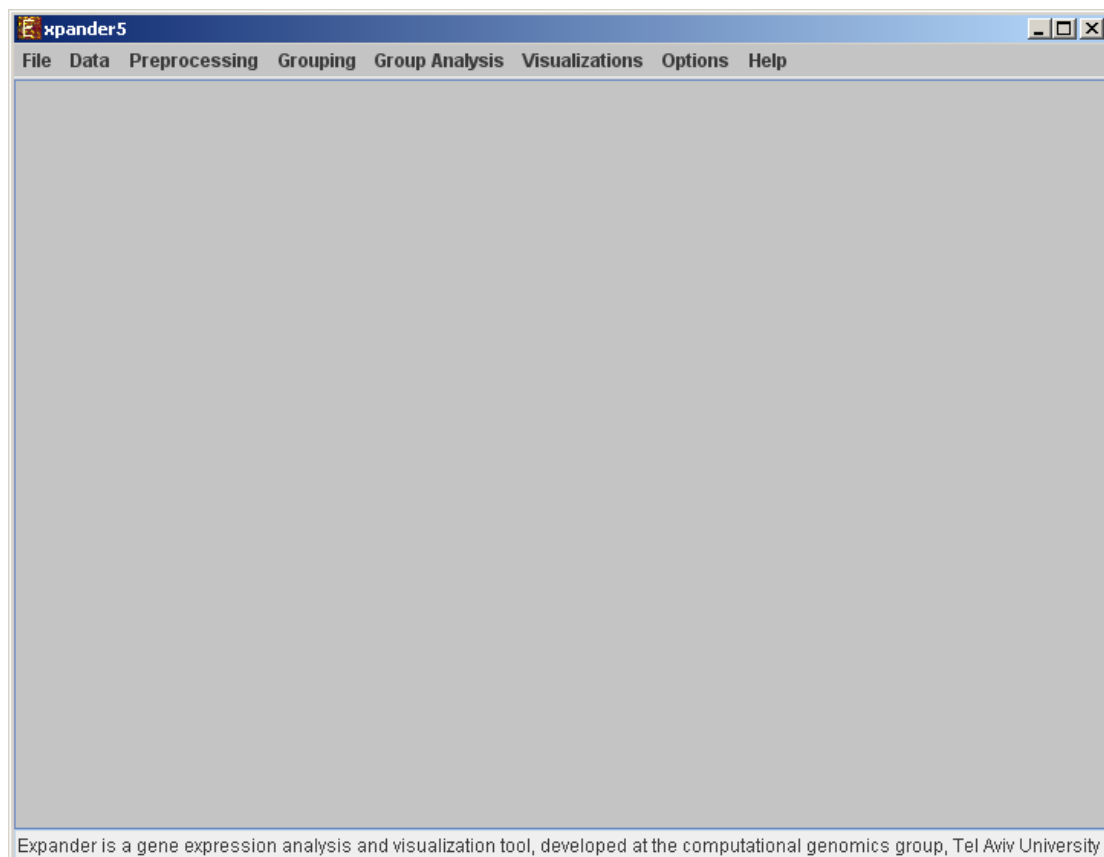

Input Data

Expander operates on two types of data:

a) **Gene expression data** – can be either relative intensities data, expected as log₂ (R/G) values data (e.g. cDNA microarrays) OR absolute intensities data, expected as positive expression levels (E.g. High-density oligonucleotide data). Oligonucleotide data can be loaded with/without detection. It can also be loaded from CEL files (If R is installed).

If one wishes to perform functional analysis or promoter analysis, an **ID conversion file** should be loaded along with the data file. The conversion file maps each probe ID (first column) in the data file into a corresponding conventional gene ID that is used in the GO annotation and TF fingerprint files that are supplied with EXPANDER. The conversion file can be loaded in the middle of the session too, by Data >> Load Conversion File.

b) **Gene sets data** – contains predefined sets of genes. In this data type, the conventional gene IDs that are used by EXPANDER in the GO annotation and TF fingerprint files are expected.

For details regarding the Gene ID convention that is used for each organism, refer to the [Supplied Files](#) section.

For details regarding the data files formats see the [File Formats](#) section.

Loading gene expression data:

Tabular Data File

To load tabular expression data, select: *File >> New Session*. From the submenu select *Expression Data >> Tabular Data File*.

When selecting *Tabular Data File*, the following dialog box will appear:

Load Tab Delimited Data

Organism: **human** Expected gene IDs: Entrez

Raw data file: **Browse**

IDs conversion file: **Browse**

Use probe IDs as gene IDs

Data type: **Relative Intensities** Data scale: **Log2**

File contains detection calls (A, M, P flags)


Set missing values to Estimate missing values with KNN

Floor value:

OK **Cancel** **Advanced**

Data type and scale are to be determined according to the input file. If the file contains missing values, these values will be estimated upon loading the data either by setting them to an arbitrary value (if the 'Set missing value to ____' option is selected) or by utilizing the KNN (K-Nearest Neighbors) method (if the 'Estimate missing values with KNN' option is selected). If the file contains Affymetrix detection calls data, the relevant check box must be checked. You may change / erase the default floor value, to which all entries that are below that value will be set (this option is available only for absolute intensities data).

Advanced Input Dialog: Upon pressing the 'Advanced' button after filling the 'Raw Data File' field, an 'Advanced Input Dialog' appears. This dialog box can be used in order to facilitate the data load of files that are not in the required format. The first few rows and columns of the data are displayed in a table, demonstrating the way the data is read by the program according to the current input values.

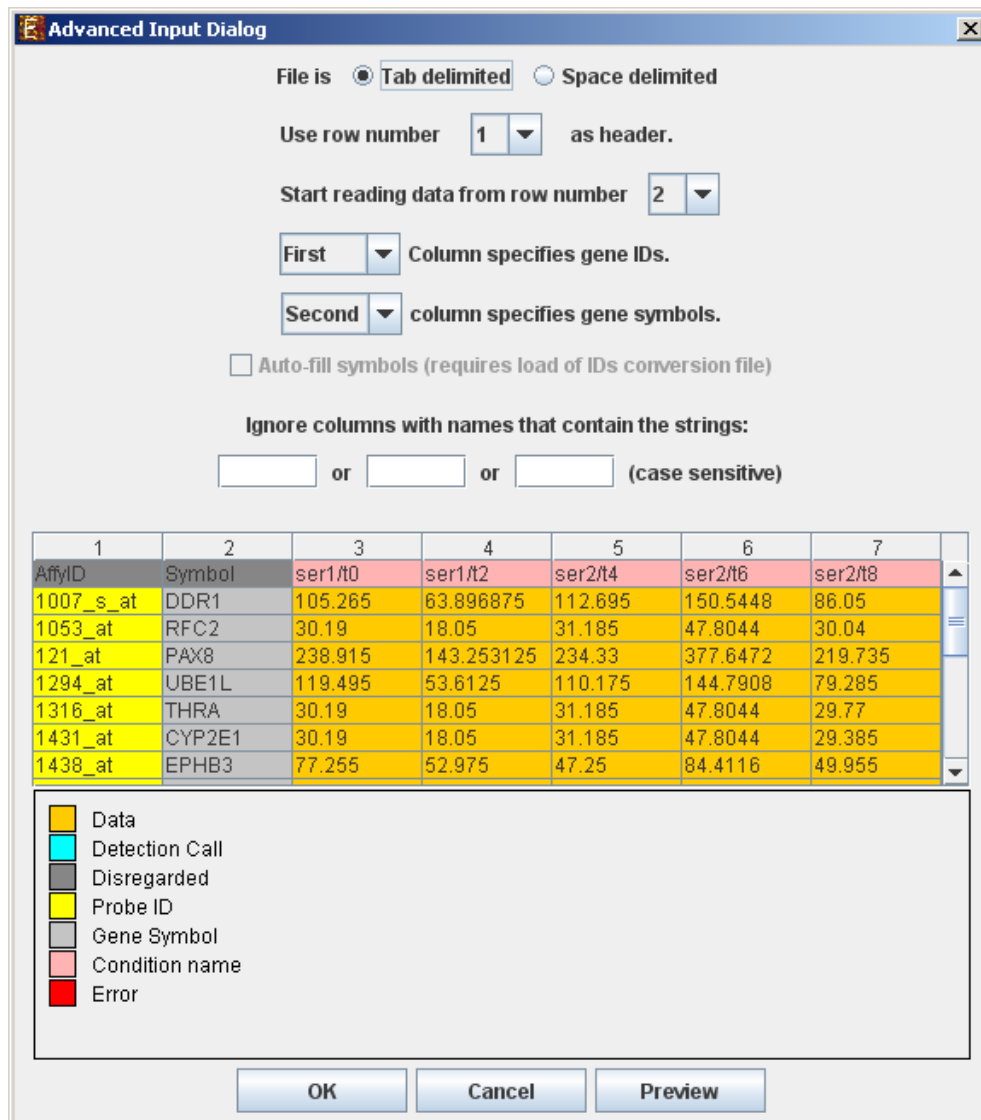

CEL Files

To load expression data from CEL files, select: *File >> New Session*. From the submenu select *Expression Data >> CEL Files*.

The load of CEL requires installation of R software (see [R External Application](#) section) along with specific packages, as detailed below. An open internet connection is also required for this operation.

Expander supports CEL files of three chip types:

1. **3' Gene Expression** - requires Bioconductor "affy" package
2. **Whole-Transcript Gene Expression** (Gene 1.0 chips) – requires "aroma.affymetrix" package.
3. **Alternative Splicing** (Exon 1.0 chips) – requires "aroma.affymetrix" package. Please note that we estimate the overall expression for the transcript, not

exon-by-exon. Therefore, this becomes 'gene data' rather than 'alternative splicing data'.

When selecting *CEL Files*, the following dialog box will appear:

The dialog box is titled "Load CEL files" and contains the following elements:

- Organism:** A dropdown menu set to "human".
- Chip type:** A dropdown menu set to "3' Gene Expression".
- Files location:** An empty text field with a "Browse" button to its right.
- Save data into file:** An empty text field with a "Browse" button to its right.
- Retrieve detection calls (A, M, P flags):** An unchecked checkbox.
- Normalization type:** Two radio buttons: "RMA" (selected) and "GC-RMA" (unselected).
- CDF:** A section containing three options:
 - Use Bioconductor default:** Selected with a radio button.
 - Import file:** An empty text field with a "Browse" button to its right.
 - Use R package:** An empty text field with a "Browse" button to its right.
- Buttons:** "OK" and "Cancel" buttons at the bottom.

Please choose the relevant organism and chip type. Then browse to the folder where the CEL files are located (*Files location*), and choose where to save the expression file resulting from the CEL files preprocessing.

Preprocessing and normalization method: The default method in Expander is RMA. However, for 3' gene expression arrays, you may select GC-RMA instead (taking into account GC-content bias). Before using GC-RMA, please make sure you have the "gcrma" R package installed (see [R External Application](#) section).

CDF environment choice: You may use the default Bioconductor CDF environment for the chips or browse to an alternative CDF package which you have already installed in R. For whole transcript and alternative splicing chips (for which there is no default Bioconductor CDF environment), you can either supply an alternative CDF package or a CDF file suitable for the chip (see links below).

Note: GC-RMA requires the probe sequence information of the chip. If you decide not to use the default Bioconductor CDF environment, and have GC-RMA as the preprocessing method, you must have the suitable probe package installed in addition to the CDF alternative package.

Links for downloading CDF files (available at [aroma.affymetrix Google group pages](#)):

- Gene 1.0 ST chips:
 - Human: [HuGene-1_0-st-v1,r3.cdf](http://hgsc.lbl.gov/cdfFiles/HuGene-1_0-st-v1,r3.cdf)
 - Mouse: [MoGene-1_0-st-v1,r3.cdf](http://hgsc.lbl.gov/cdfFiles/MoGene-1_0-st-v1,r3.cdf)
 - Rat: [RaGene-1_0-st-v1,r3.cdf](http://hgsc.lbl.gov/cdfFiles/RaGene-1_0-st-v1,r3.cdf)
- Exon 1.0 ST chips (we recommend to download a recent 'core' cdf file):
 - Human: http://bcgc.lbl.gov/cdfFiles/HuEx-1_0-st-v2,A20071112,EP/HuEx-1_0-st-v2,coreR3,A20071112,EP.cdf
 - Mouse: http://bioinf.wehi.edu.au/folders/mrobinson/CDF/MoEx-1_0-st-v1,coreR1,A20080718,MR.cdf
 - Rat: http://bioinf.wehi.edu.au/folders/mrobinson/CDF/RaEx-1_0-st-v1,coreR1,A20080718,MR.cdf


Note that new CDF files may become available before we get to adjust our links accordingly.

For Exon chips, there are other CDF files available from aroma.affymetrix group. The links we supplied are for core CDF files, which allow faster analysis than the 'extended' or 'full' CDF files.

If Expander cannot find your R software, a window will appear, asking you to specify its location. Please browse to the location of your R software. In Windows, R.exe file is likely to be located in the 'bin' folder of R software. In Linux, you may type 'which R' in the command line to find R path. If you have a few versions of R installed, please make sure to point Expander to a version in which the Bioconductor "affy" and/or "aroma.affymetrix" packages have been installed.

Once the CEL files preprocessing is done, a corresponding tabular data file is generated and a 'Load Study' dialog will appear, as in loading [Tabular Data](#).

After loading a gene expression data set, a 'Session Data' display tab is added to the main window (see example below). It contains information regarding the raw data file, a box plot chart, and an expression matrix visualization of the raw data. If detection calls exist in the data file, their statistics for each probe appear in 3 columns in the heat maps (expression matrices), in a scale between 0 and 1, corresponding to the relative part of each of the detection calls (P, M and A). The detection calls statistics for each condition are displayed in a separate tab in two tables (one for the raw data and another for the preprocessed data) and are presented in percent.

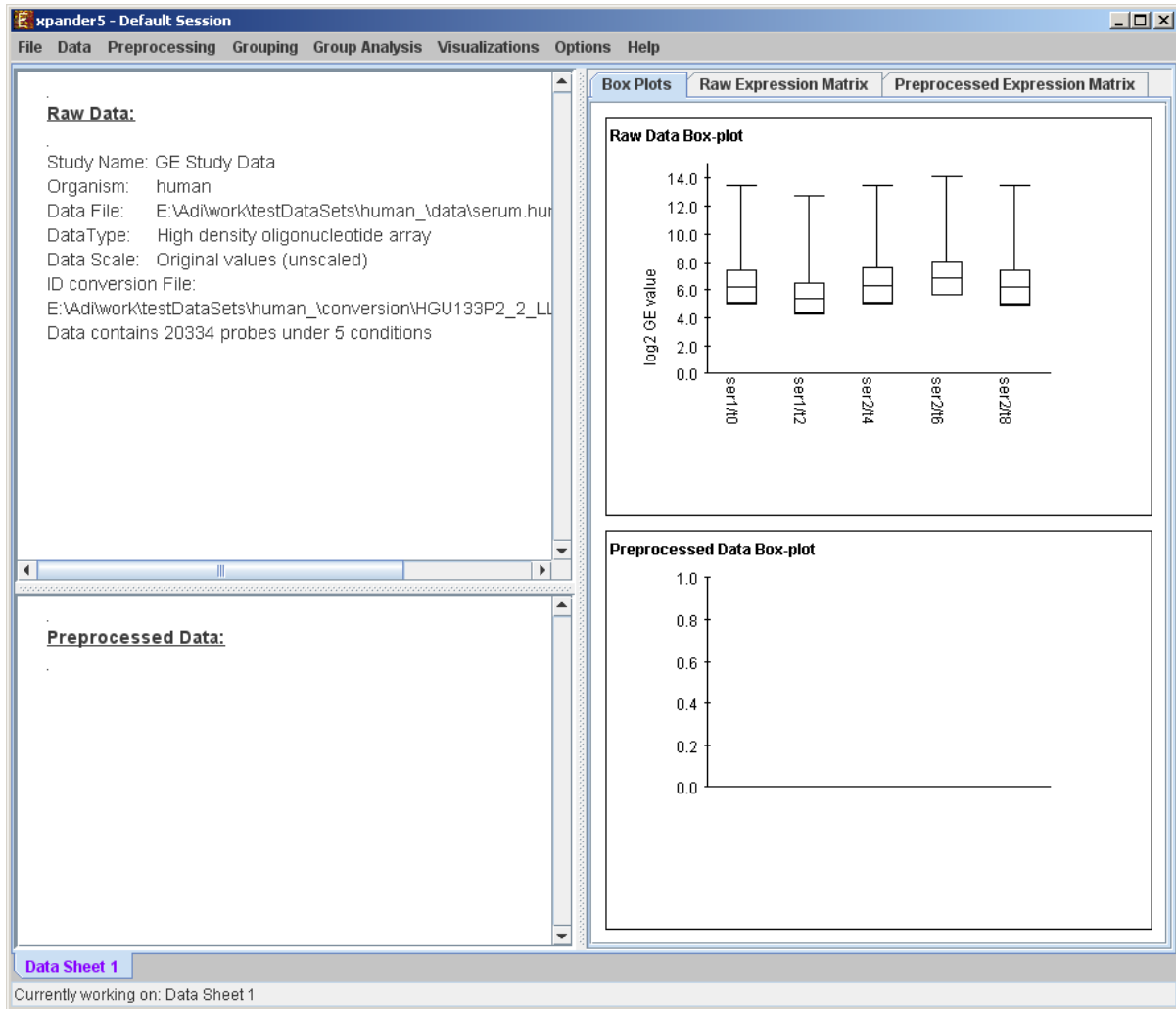

Loading gene sets:

Gene Sets:

To load Gene Sets data select *File>>New Session*. From the submenu select *Gene Sets*.

The following dialog box will appear:

The 'Load gene sets' dialog box contains the following fields and options:

- Organism:** A dropdown menu set to 'human'.
- Expected Gene IDs:** A dropdown menu set to 'Entrez'.
- Data name:** A text input field containing 'Gene Sets Data'.
- File Name:** A text input field with a dropdown arrow and a 'Browse' button next to it.
- Delimiters:** Two radio buttons: 'Tab delimited' (selected) and 'Space delimited'.
- Buttons:** 'OK' and 'Cancel' buttons at the bottom.

For details regarding the data files formats see the [File Formats](#) section.

After loading a gene sets data, a 'Session Data' display tab is added to the main window (see example below). It contains information regarding the data file, and a table describing the different sets (set number, name and size). Set names can be modified, by editing the corresponding cell in the table. Upon clicking on a row in the table, the corresponding set pane appears on the right. It contains a list of the genes in the set and a view of their chromosomal locations. If a network file has been loaded (via Data>>Load Network), the sub-graph, induced by the set is displayed as well.

The screenshot shows the xpander4 software interface. The main window is titled 'xpander4 - Default Session'. The menu bar includes File, Data, Preprocessing, Grouping, Group Analysis, Visualizations, Options, and Help. The main area is divided into two panes. The left pane, titled 'Gene sets data:', shows the following information:

Name: Gene Sets Data
 Organism: mouse
 Gene sets loaded from file: E:\Adi\work\testDataSets\mouse\ttt_rmm_CWM.txt

Sets table:

ID	Name	Size
13	13	1581
18	18	1389
43	43	1386
14	14	1322
4	4	1318
6	6	1213
49	49	1198
35	35	1189
41	41	1167
20	20	1157
21	21	1086
15	15	1057
9	9	1047
45	45	1040
50	50	1032
30	30	993
28	28	985
27	27	975
24	24	960
46	46	918
1	1	888
47	47	888
8	8	885
19	19	807
40	40	727

The right pane, titled '41 Set content', shows a list of genes with their Gene ID and Gene Symbol:

Gene ID	Gene Symbol
74088	0610012H03Rik
20707	Serpinc9c
268905	Krtap13-1
207798	Gramd1c
71743	Coasy
16667	Krt17
69464	2300006N05Rik
68371	Pblid
16821	Lcn4
18128	Notch1
215928	BC021785
232493	Gys2
30805	Slc22a4
26378	Decr2
14408	Gabrr1
75564	1700027N10Rik
230576	Ttc22
59083	Fatub
14473	Gc
235505	Cd109
68352	0610012D14Rik
114142	Foxp2
71775	1300017J02Rik
22157	Tulp1
75578	2310009E04Rik
74188	Prlnc3
18010	Neu1
20720	Serpine2
212943	BC023892
76585	2310069N01Rik
70445	Cd248
14725	Lrp2
58875	Hibadh
14865	Gstm4
56316	Ggcx
433180	EG433180
12499	Entpd5
68239	Krt42
76500	Ihpk2
228802	BC018465
231532	Arhaap24
16664	Krt14

At the bottom of the window, a status bar reads: 'Expander is a gene expression analysis and visualization tool, developed at the computational genomics group, Tel Aviv University'.

Preprocessing GE Data

The following preprocessing schemes can be performed using EXPANDER:

- 1) **Flooring:** setting all expression values that are below a certain threshold (set by the user) into that threshold. This can be performed through the 'Data Load' dialog box, and is available only for absolute intensities data.
- 2) **Merging conditions** (*Preprocessing >> Merge conditions*): merging a selected set of condition profiles (columns) in the dataset into one profile, in which each entry holds the average value of the merged entries.
- 3) **Merging probes according to gene ID** (*Preprocessing >> Merge Probes by Gene ID*): automatically shrinks the matrix so that all rows of probes from the same gene are merged into one average row, identified by the corresponding gene ID.
- 4) **Normalization:** required in order to remove systematic variation, i.e. variation arising from reasons other than biological differences between RNA samples.

Expander performs normalization only for absolute intensities data, since it is assumed that the relative intensities data (e.g. cDNA microarrays) is already normalized, as it is input after performing log ratio ($\log_2 R/G$).

Normalization can be performed using the following schemes:

- a) **Quantile normalization** (*Preprocessing >> Normalization >> Quantile*), in which the whole data is used.
- b) **Non-linear baseline normalization** (*Preprocessing >> Normalization >> Non Linear Baseline*), which uses a baseline array (can be selected by the user). In this scheme a normalization function is calculated using pseudo Loess regression of the M vs. A scatter plot. The subset of genes that are used to evaluate the normalization function can be set to 'all genes' (recommended when most genes in the dataset are expected to be constantly expressed) or a 'rank invariant set' of genes (recommended when there can be a large number of differentially expressed genes).

For more details regarding the normalization schemes see the [References](#) section.

- 5) **Condition filtration:** the conditions used in the analysis can be manually filtered by selecting: *Preprocessing >> Filter Conditions*. This will bring up a dialog box in which the user can select the required conditions from a list.
- 6) **Gene (probe) filtration:** can be performed in order to filter out some of the constantly expressed genes, and perform downstream analysis on a smaller informative subset of the genes.

Probe filtration can be performed using the following schemes:

- a) **t-Test** (*Preprocessing >> Filter Probes >> t-Test*): When using this method, only probes that demonstrate differential expression between two condition subsets are selected.
- b) **SAM - Significance Analysis of Microarray** (*Preprocessing >> Filter Probes >> SAM*): selects probes that demonstrate differential expression between conditions subsets. You may choose 2 or more subsets (multi-class tests are supported). This method uses permutations to get an 'empirical' estimate for the FDR of the reported differential genes (for details see the [References](#) section). Before using SAM, please make sure you have **R software along with the "samr" package** installed (see [R External Application](#) section).
- c) **Fold Change** (*Preprocessing >> Filter Probes >> Fold Change*): when using this method only genes that are over/under expressed by at least n fold in at least k arrays are selected (n and k are determined by the user). The fold change can be calculated in relation to (a) a selected baseline array (b) the minimal expression value of the gene OR (c) the reference value when working on relative intensities (depending on the user's selection).
- d) **Variation** (*Preprocessing >> Filter Probes >> Variation*): In this method, the k most variant genes are selected (k is determined by the user). Variance is used to measure variation for relative intensities data, and Coefficient of Variation is used to measure variation for absolute intensities data.
- e) **Detection calls** (*Preprocessing >> Filter Probes >> Detection calls*): in this method probes/genes are filtered according to the number of expression signals for which the detection call is 'P' (Present). It can only be operated if the data file contains detection info.

- f) **Load Probe Subset** (*Preprocessing >> Filter Probes >> Load Probe Subset*): the filtered set is loaded from an external txt file (for details regarding the format please see the [File Formats](#) section).
- 7) **Standardization**: When expression values between different genes are very different, but general expression patterns are similar (high Pearson Correlation values), we would expect to see this similarity when looking on a pattern display. Since the absolute values of expression are different, a manipulation is required, in order to view the patterns on the same scale. This manipulation is called standardization.

Standardization can be performed using the following schemes:

- a) **Mean 0 and Variance 1** (*Preprocessing >> Standardization >> Mean 0 and Variance 1*) – normalizes each expression pattern to have a mean of 0 and a variance of 1. This method is appropriate in most cases when working on genes.
- b) **Log data** (*Preprocessing >> Standardization >> Log Data*) – Performs log₂ operation on each entry.
- c) **Fixed norm** (*Preprocessing >> Standardization >> Fixed Norm*) - normalizes each expression pattern to have a fixed norm i.e. expression levels are divided by the norm of that expression vector (the root of sum of squares of that vector). This method is appropriate when different mean values or variances are expected for different patterns (e.g. when working on conditions and expecting larger variance in later phases of a response).
- d) **Divide by Base** (*Preprocessing >> Standardization >> Divide by Base*) – Divides each entry in a profile (a column) by the corresponding entry in the profile of a selected base condition. This can be done for all conditions or for subsets of the conditions.

After performing a preprocessing operation, the information regarding the operation is added to the 'Preprocessed Data' section in the 'Session Data' tab. In addition, the 'Preprocessed Data box plot' and 'Preprocessed Expression Matrix' are automatically updated according to the new values in the data.

File **Preprocessing** **Grouping** **Group Analysis** **Visualizations** **Options** **Help**

Raw Data:

Study Name: Data Sheet
 Organism: human
 Data File: D:\adi_work\testDataSets\human\data\serum.human.txt
 Data Type: High density oligonucleotide array
 Pubmed ID: None
 ID conversion File:
 D:\adi_work\testDataSets\human\conversion\HGU133P2_2_LLids.txt
 Data contains 20334 probes under 5 conditions

Preprocessed Data:

1) Quantile Normalization
 2) Fold Change Filter:
 fold change: 3.0
 entries exp: 1
 Base condition: t0
 survived filter: 750
 3) Standardization:
 Type: mean 0 and std 1

Preprocessed Expression Matrix


Box Plots **Expression Matrix**

Raw Data Box-plot

Preprocessed Data Box-plot


Data Sheet 1

Currently working on: Data Sheet 1

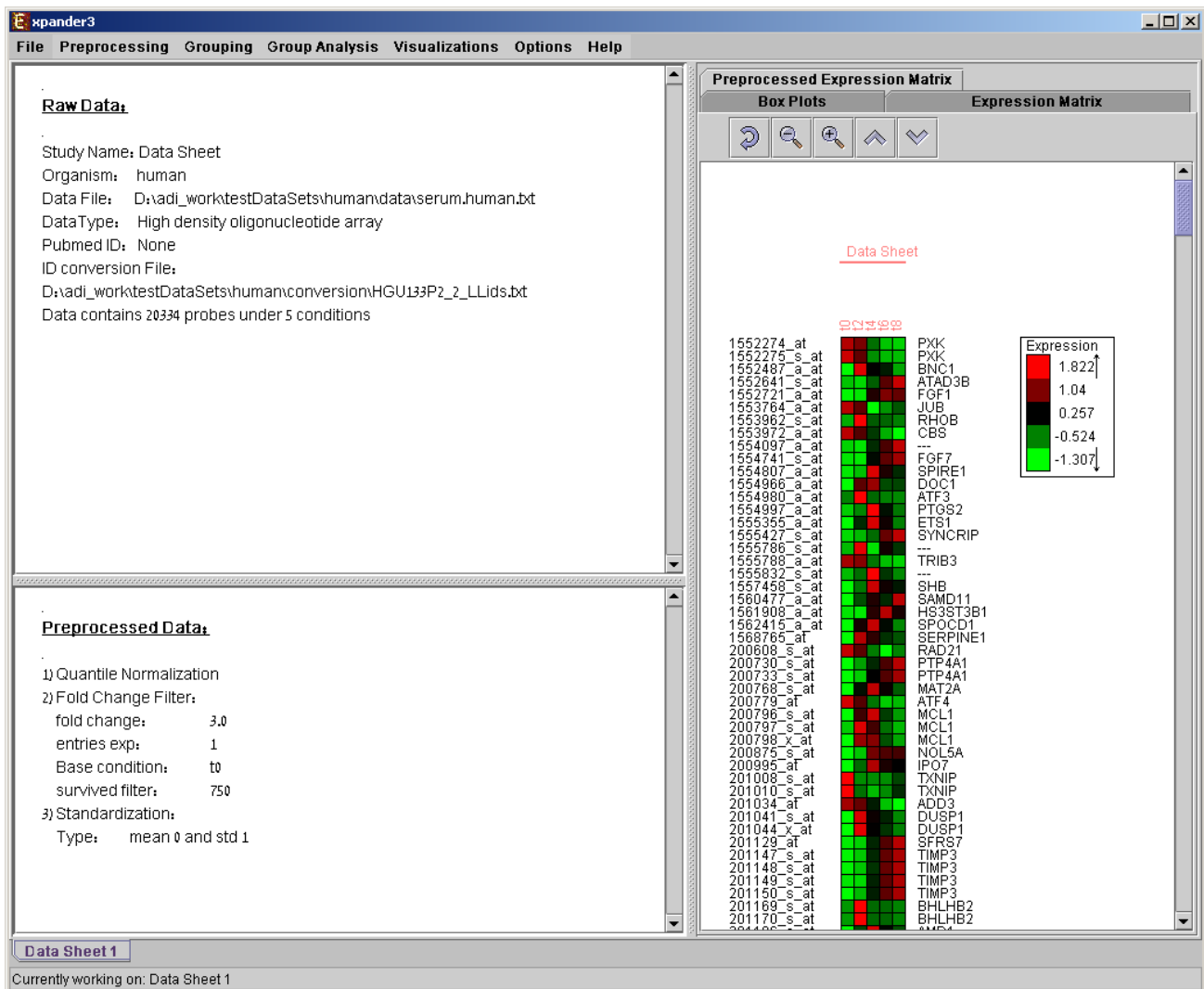

Upon selecting *Preprocessing* >> *Undo* the data is changed to be as it was before the most recent preprocessing operation was performed, and the corresponding information is removed from the 'Preprocessed Data' section. The 'Preprocessed Data box plot' and 'Preprocessed Expression Matrix' are automatically updated accordingly.

All the above operations can be performed before running further analysis on the data and generating displays. When attempting to perform further preprocessing operations after analysis results and visualizations have been generated, the following dialog box appears:

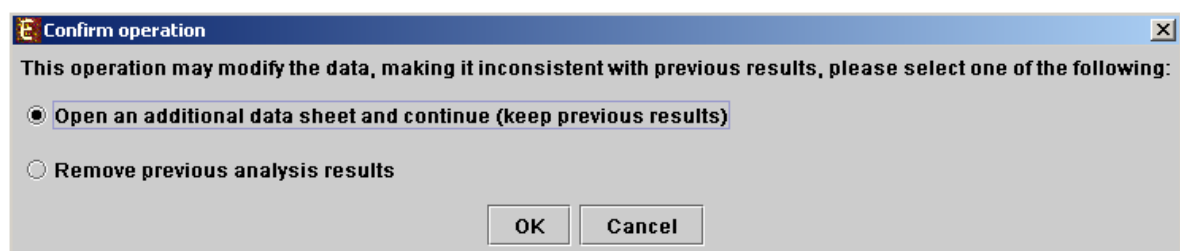


Upon choosing to open an additional data sheet, a new data set view tab called 'Data Sheet 2' is added to the main frame. The title of this tab is highlighted (colored in purple), indicating that it is now the active data sheet (i.e. all further operations refer to this data sheet). The active data sheet is automatically changed according to the selected (front) visualization tab.

Preprocessed gene expression data can be saved to a file at any time by selecting *Preprocessing >> Save Preprocessed Data*. The data is written in the same format defined for input GE data.


Viewing Data Plots

Expander provides two types of scatter plots visualizations that can be operated via *Preprocessing >> Normalization >> View Scatter Plots*.

Simple plot - Displays a scatter plot of two arrays (selected by the user), in which the i th point (x_i, y_i) represents the expression value (log expression for un-logged data) of the i^{th} gene in one array vs. the other. For normalized data, points should be located around the $y=x$ line (marked on the scatter plot).

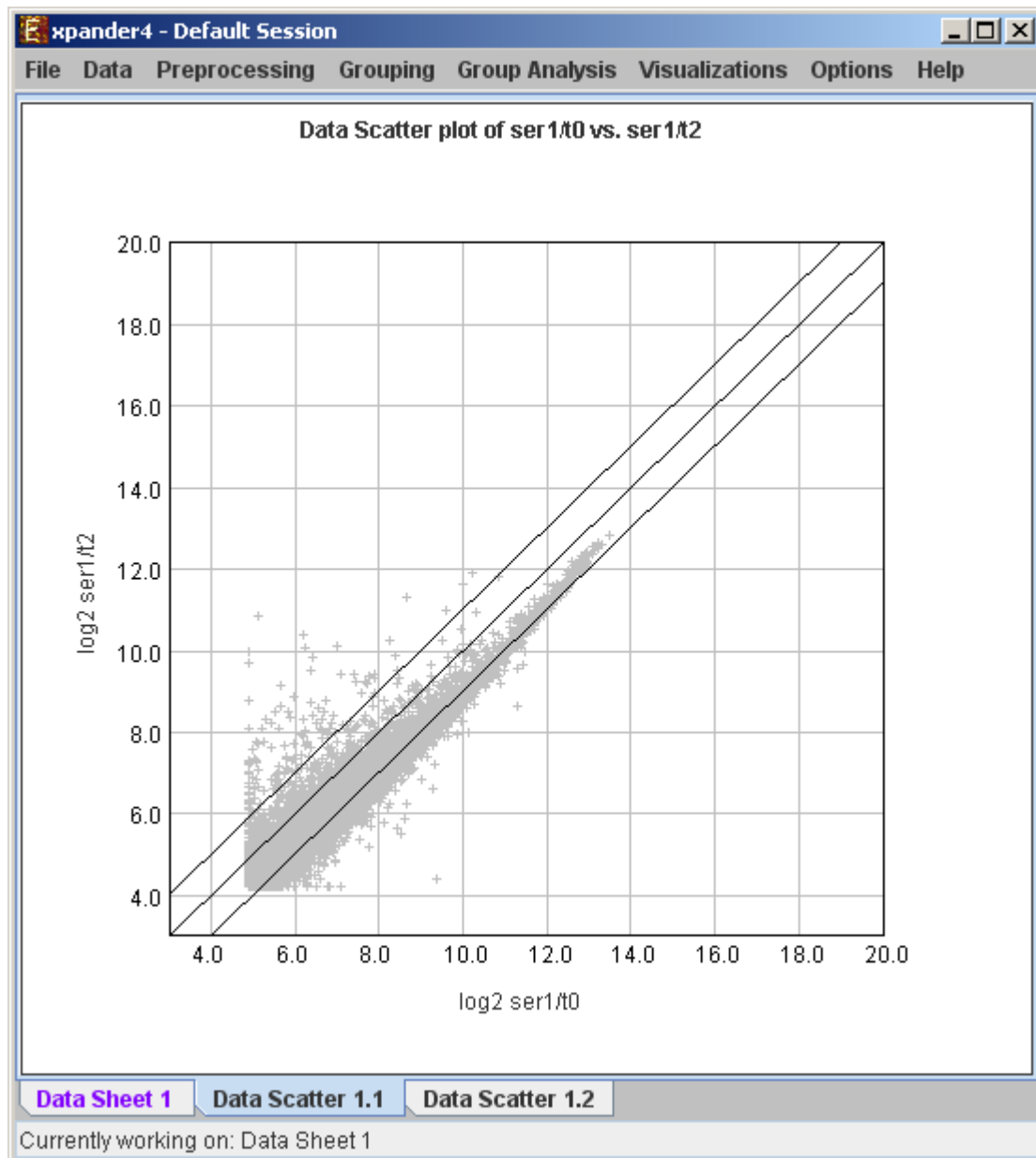

M vs. A plot (available only for absolute intensities data) - Displays a scatter plot in which each point (A_i, M_i) represents the log intensity difference of the i th probe in the two arrays (selected by the user) vs. the average log value of these intensities.

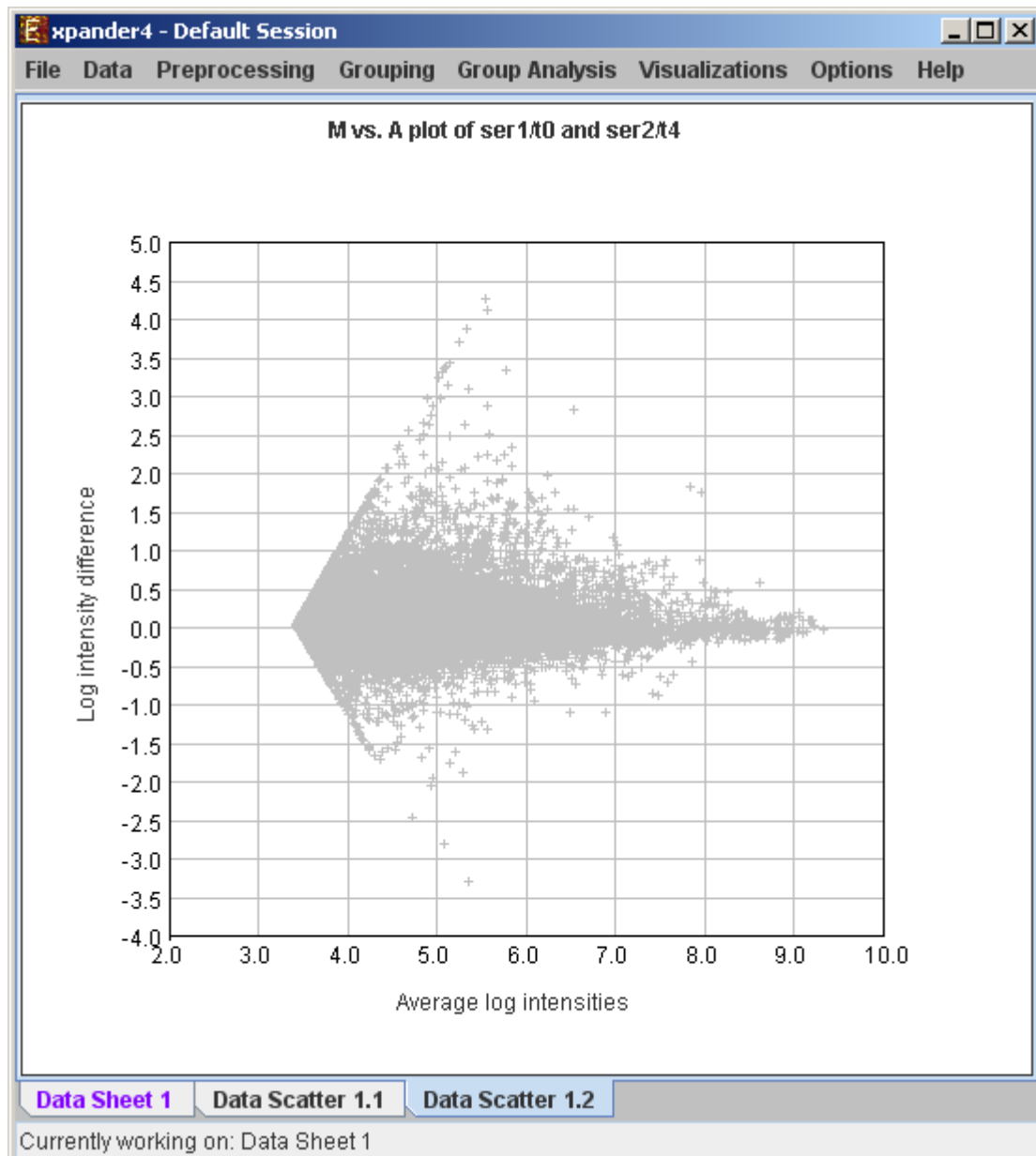

Clustering GE Data

The goal of clustering is to partition the genes into distinct sets such that genes that are assigned to the same cluster should have similar expression patterns, while genes assigned to different clusters should have non-similar expression patterns. Usually there is no one solution that is the 'true' mathematical solution for this problem, but a good clustering solution should have two merits:

- (1) High homogeneity (average similarity between genes from the same cluster).
- (2) High separation (average distance/dissimilarity between genes from different clusters).

After operating one of the clustering algorithms a clustering results view appears. The view contains information about the solution and its quality including the method and parameters that were used to obtain it, number of clusters, number of singletons (probes that were not assigned to any cluster), overall homogeneity and separation, as well as the size and homogeneity of each cluster. This summary can be used to compare different solutions.

In order to apply a clustering algorithm to the data, select the required algorithm from the *Grouping >> Clustering* menu (options are: **KMeans**, **CLICK**, **SOM**). You can also load an existing clustering solution from a file by selecting the **Load Solution** option from the *Clustering* menu (For details regarding the clustering solution file format, refer to the [File Formats](#) section). For details about **Hierarchical clustering**, please go to [Hierarchical Clustering](#).

The **CLICK** algorithm is not designed to find clusters under the size of 15 probes, so it might fail in clustering small datasets.

Fill the required input data in the algorithm input dialog box and press the 'Ok' button. The parameters required for each method are as follows:

Algorithm	Required parameters
KMeans	Expected number of clusters.
SOM	Grid width, grid length (width*length >= number of clusters) and number of iterations.
CLICK	Homogeneity value (0-1): allows the user control over the homogeneity of the resulting clustering, i.e. the average similarity between elements in the same cluster. This parameter serves as a threshold in various steps in the algorithm, including the definition of cluster kernels, singleton adoptions and kernel merging. The default value for this parameter is the estimated homogeneity of the true clustering. The higher the value assigned to this parameter the tighter the resulting clusters.

Details about the algorithms can be obtained through the relevant articles in the [References](#) section.

After clustering is performed, a clustering solution visualization tab is added to the main window. It contains the following views:

Information regarding the clustering algorithm, number of clusters, number of un-clustered elements (singletons), and numerical measures of the clustering quality, including:


- a) Overall average homogeneity - calculated as the average value of similarity between each element and the center of the cluster to which it has been assigned, weighted according to the size of the cluster.
- b) Overall average separation – calculated as the average similarity between mean patterns of different clusters, weighted according to their sizes.


- c) Clusters table - contains the number, name (label), size and homogeneity of each cluster. The name of a cluster can be changed by editing the corresponding cell in the table.

Mean Patterns of all clusters with error bars (± 1 STD).

A table listing all condition titles and their corresponding number used in the patterns display. Upon selecting a row in this table, the corresponding column in each of the mean pattern plots is marked.

Upon selecting a cluster (from the clusters table or from the mean patterns view), the corresponding cluster pane is displayed on the right. It contains a list of probes, probe patterns, expression matrix (heat map) and the chromosomal locations of the genes. Similarity matrices for probes within the cluster as well as for conditions are also displayed in this tab, if the relevant options in the display settings are selected (see the [Settings](#) section). If a network file has been loaded (via Data>>Load Network), the sub-graph, induced by the cluster is also displayed in the cluster pane.

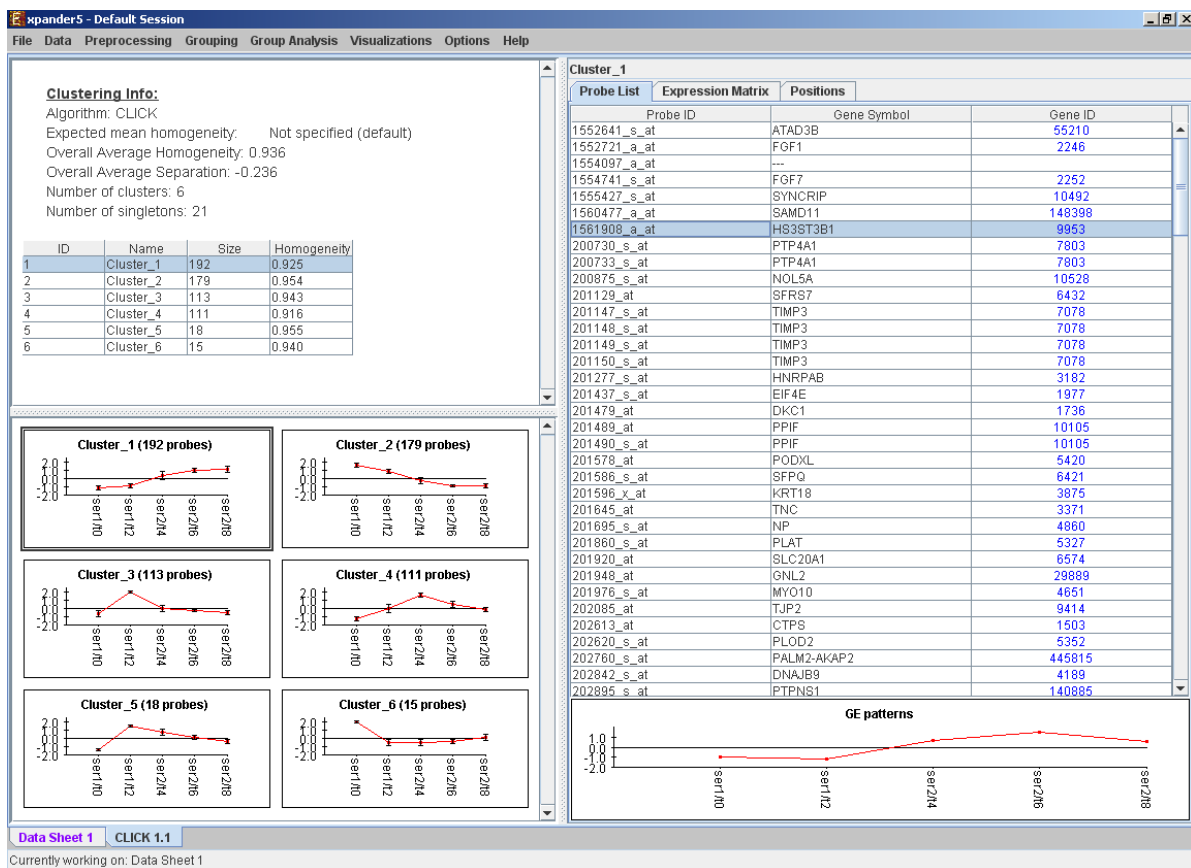

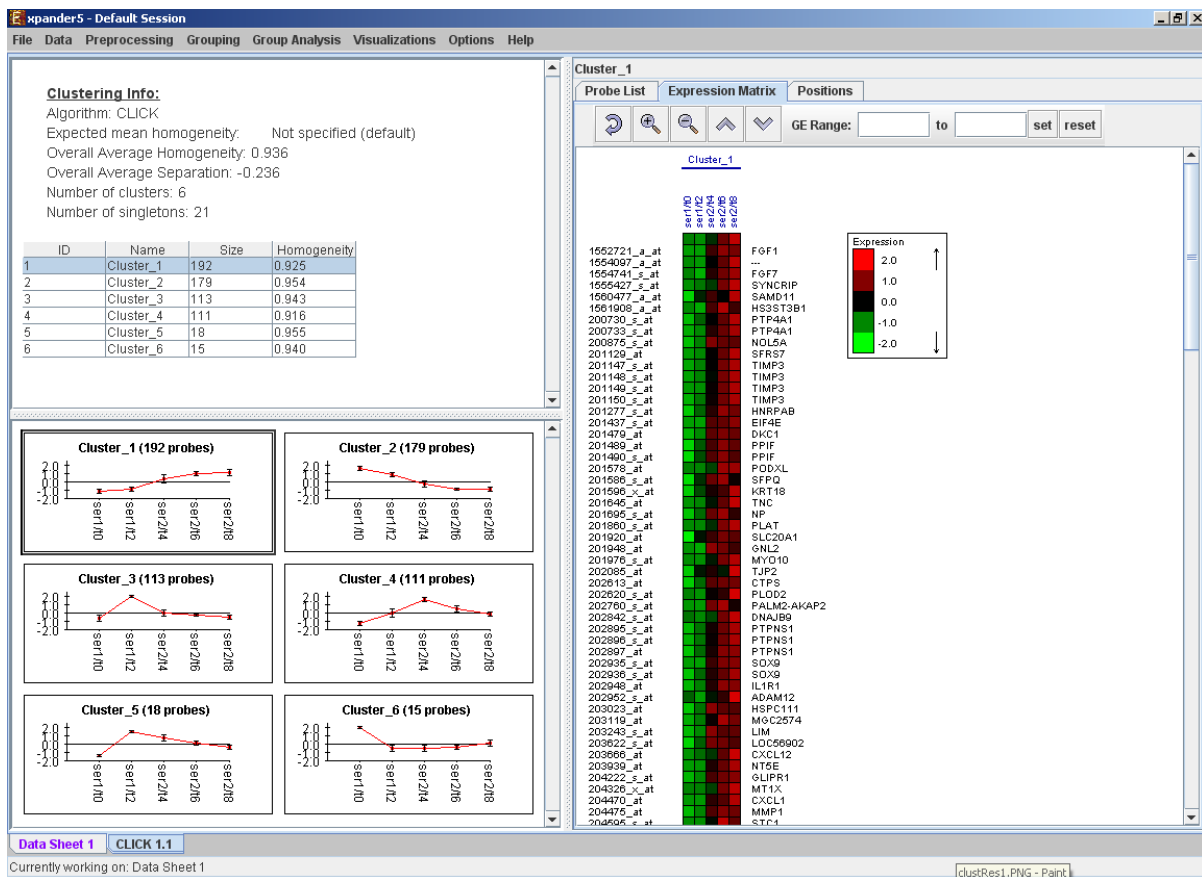

After performing group analysis (for details see the [Group Analysis Tools](#)), if enrichment has been detected in the selected cluster, the corresponding histogram and analysis information are added to the single cluster view.

A clustering solution can be saved using the *Grouping >> Clustering >> Save Solution*, and reloaded using the *Grouping >> Clustering >> Load Solution*.

Defining a group according to similarity to a selected probe


This can be done by selecting *Grouping >> Group by Pattern Similarity*. An input dialog box allows setting the similarity measure (Pearson correlation, Spearman correlation or Euclidean distance) and reference probe ID as well as the expected group size. The visualization for this operation is similar to the clustering results visualization (described above).

Biclustering GE Data

Biclustering is clustering of both genes and conditions of the data into subgroups that are not necessarily disjoint. It enables the user to detect genes that are co-regulated in only a subgroup of the conditions, and does not force genes to belong exclusively to one cluster. It is useful when working on datasets which contain a large number of conditions.

Biclustering is performed by Expander using the SAMBA algorithm (for details see the [References](#) section).

In order to apply the SAMBA biclustering algorithm to the data select *Grouping>>Bi-Clustering>>SAMBA*. The following dialog box will appear:

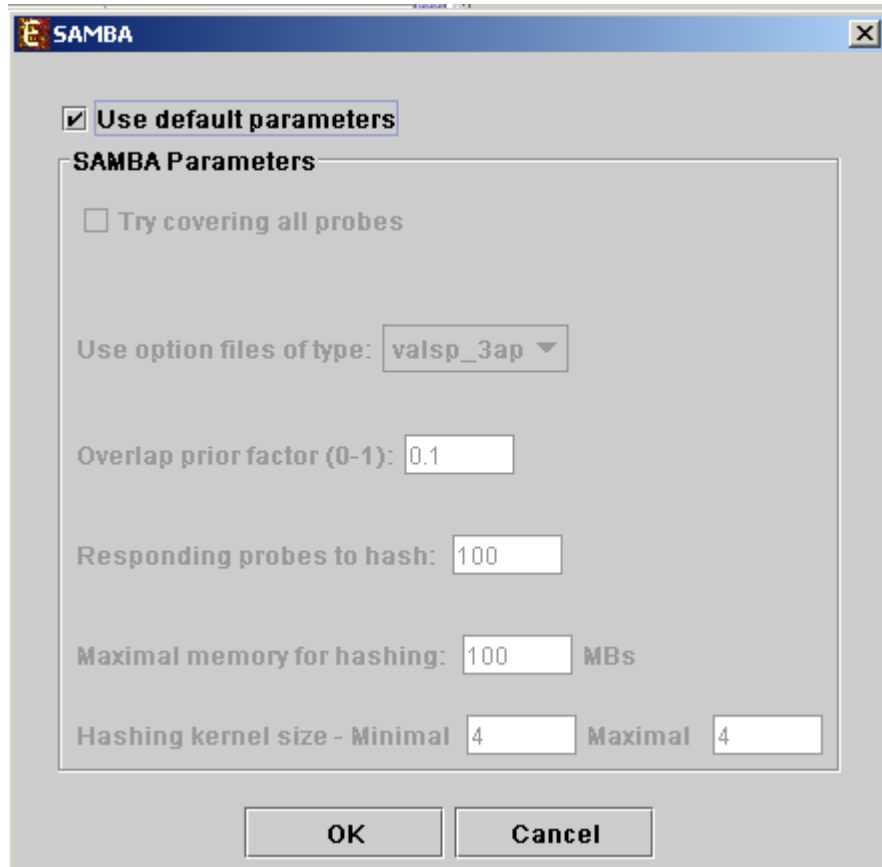

It enables the configuration of some of the parameters for the algorithm. The following table specifies the different parameters that can be set via this dialog box:


Field	Description																				
Use default parameters	When checked, biclustering parameters (described below) are set automatically (this option is recommended unless the user is familiar with the parameters).																				
Option files type	The user can select one out of 6 options. The following table describes the advantages and disadvantages of each option: <table border="1" style="margin-left: 20px;"> <thead> <tr> <th>Option name</th> <th>fast performance</th> <th>less memory required:</th> <th>Flexible</th> <th>Robust- can handle normalization problems and non gene-expression data</th> </tr> </thead> <tbody> <tr> <td>valsp_1</td> <td style="text-align: center;">+</td> <td style="text-align: center;">+</td> <td style="text-align: center;">-</td> <td style="text-align: center;">-</td> </tr> <tr> <td>valsp_2</td> <td style="text-align: center;">0</td> <td style="text-align: center;">0</td> <td style="text-align: center;">0</td> <td style="text-align: center;">-</td> </tr> <tr> <td>valsp_3</td> <td style="text-align: center;">-</td> <td style="text-align: center;">-</td> <td style="text-align: center;">+</td> <td style="text-align: center;">-</td> </tr> </tbody> </table>	Option name	fast performance	less memory required:	Flexible	Robust- can handle normalization problems and non gene-expression data	valsp_1	+	+	-	-	valsp_2	0	0	0	-	valsp_3	-	-	+	-
Option name	fast performance	less memory required:	Flexible	Robust- can handle normalization problems and non gene-expression data																	
valsp_1	+	+	-	-																	
valsp_2	0	0	0	-																	
valsp_3	-	-	+	-																	

	valsp_1ap	+	+	-	+
	valsp_2ap	0	0	0	+
	valsp_3ap	-	-	+	+
	We recommend the valsp_3ap option (set as default), since it is very flexible, and produces good results also for data that was not normalized properly or for non gene-expression data.				
Always cover all genes	When checked, the solution will cover each gene at least once (each gene will be included in one or more biclusters).				
Always cover all conditions	When checked, the solution will cover each condition at least once (each condition will be included in one or more biclusters). Unchecking this option will cause a reduction in the number of biclusters, and the algorithm will run faster.				
Overlap prior factor	Can take values between 0 and 1, describes extent of overlap that is permitted between two different biclusters in the same solution. The higher this parameter is, the more strict the algorithm will be regarding adding a new bicluster (will require less overlap between the new bicluster and the existing ones).				
Number of responding genes to hash	Can take values between 1 and the number of genes in the dataset. Default value is set to 100 (recommended unless data set size < 100). Has impact over the hashing stage in the algorithm.				
Maximum hash size (in MB)	Described the maximum memory size that can be used for the hashing part of the algorithm (the whole algorithm will take up about twice this size of memory).				
Maximum hash size	This parameter determines the number of condition kernel options that are tested and scored in the hashing stage. It can take values from 1 to 7. The default value is 4. In datasets with many conditions raising this number will significantly increase the algorithm run time (may also produce better results).				
Minimum hash size	This parameter determines the minimal size of condition kernel in the hashing stage. It can take values from 1 to 7 and must be <= Maximum hash size. The default value is 4.				


Upon clicking 'OK' in the dialog box, the SAMBA algorithm is operated on the dataset.

After biclustering is performed a biclustering solution visualization tab is added to the main window. It contains the following views:

- a) Information regarding the biclustering algorithm, and number of resulting biclusters.
- d) Biclusters table – contains the following information for each bicluster: serial number, name, score, number of probes genes and number of conditions. The name of a bicluster can be changed by editing the corresponding cell in the table. The score is given by the SAMBA algorithm and is size-dependent, thus, it is not recommended to use it to compare the quality of two biclusters of different sizes. The table can be filtered to display a subset of the biclusters

by clicking on the 'Filter' () button in the toolbar. Filtering can be performed according to: Score, number of probes and number of conditions.

Upon selecting a bicluster (from the biclusters table), the corresponding pane is displayed on the right. It contains a list of probes, probe patterns, expression matrix (heat map) and the chromosomal locations of the genes. Similarity matrices for probes within the cluster as well as for conditions are also displayed in this tab, if the relevant options in the display settings are selected (see the [Settings](#) section). If a network file has been loaded (via Data>>Load Network), the sub-graph, induced by the cluster is also displayed in the cluster pane.

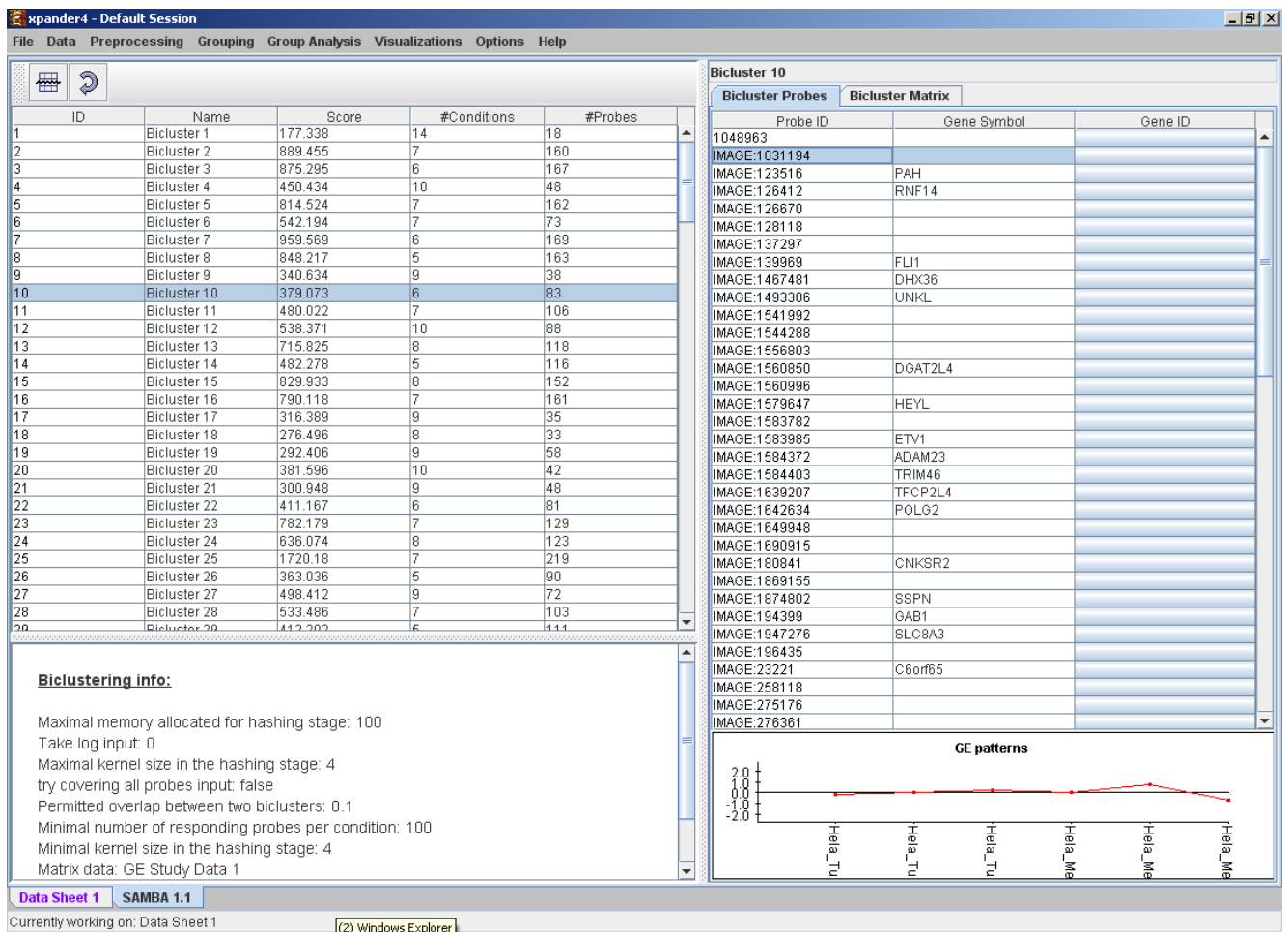

Expander4 - Default Session
File Data Preprocessing Grouping Group Analysis Visualizations Options Help

ID	Name	Score	#Conditions	#Probes
1	Bicluster 1	177.338	14	18
2	Bicluster 2	889.455	7	160
3	Bicluster 3	875.295	6	167
4	Bicluster 4	450.434	10	48
5	Bicluster 5	814.524	7	162
6	Bicluster 6	542.194	7	73
7	Bicluster 7	959.569	6	169
8	Bicluster 8	848.217	5	163
9	Bicluster 9	340.634	9	38
10	Bicluster 10	379.073	6	83
11	Bicluster 11	480.022	7	106
12	Bicluster 12	538.371	10	88
13	Bicluster 13	715.825	8	118
14	Bicluster 14	482.278	5	116
15	Bicluster 15	829.933	8	152
16	Bicluster 16	790.118	7	161
17	Bicluster 17	316.389	9	35
18	Bicluster 18	276.496	8	33
19	Bicluster 19	292.406	9	58
20	Bicluster 20	381.596	10	42
21	Bicluster 21	300.948	9	48
22	Bicluster 22	411.167	6	81
23	Bicluster 23	782.179	7	129
24	Bicluster 24	636.074	8	123
25	Bicluster 25	1720.18	7	219
26	Bicluster 26	363.036	5	90
27	Bicluster 27	498.412	9	72
28	Bicluster 28	533.486	7	103
29	Bicluster 29	412.202	6	111


Bicluster 10

Probe ID	Gene Symbol	Gene ID
1048963		
IMAGE:1031194		
IMAGE:123516	PAH	
IMAGE:126412	RNF14	
IMAGE:126670		
IMAGE:128118		
IMAGE:137297		
IMAGE:139969	FLI1	
IMAGE:1467481	DHX36	
IMAGE:1493306	UNKL	
IMAGE:1541992		
IMAGE:1544288		
IMAGE:1556803		
IMAGE:1560850	DGAT2L4	
IMAGE:1560996		
IMAGE:1579647	HEYL	
IMAGE:1583782		
IMAGE:1583985	ETV1	
IMAGE:1584372	ADAM23	
IMAGE:1584403	TRIM46	
IMAGE:1639207	TFCP2L4	
IMAGE:1642634	POLG2	
IMAGE:1649948		
IMAGE:1690915		
IMAGE:180841	CNKSR2	
IMAGE:1869155		
IMAGE:1874802	SSPN	
IMAGE:194399	GAB1	
IMAGE:1947276	SLC8A3	
IMAGE:196435		
IMAGE:23221	C6orf85	
IMAGE:258118		
IMAGE:275176		
IMAGE:276361		



Biclustering info:

- Maximal memory allocated for hashing stage: 100
- Take log input: 0
- Maximal kernel size in the hashing stage: 4
- try covering all probes input: false
- Permitted overlap between two biclusters: 0.1
- Minimal number of responding probes per condition: 100
- Minimal kernel size in the hashing stage: 4
- Matrix data: GE Study Data 1

GE patterns

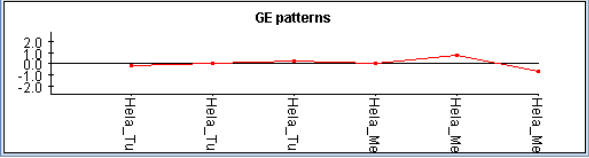

Data Sheet 1 SAMBA 1.1
Currently working on: Data Sheet 1 (2) Windows Explorer

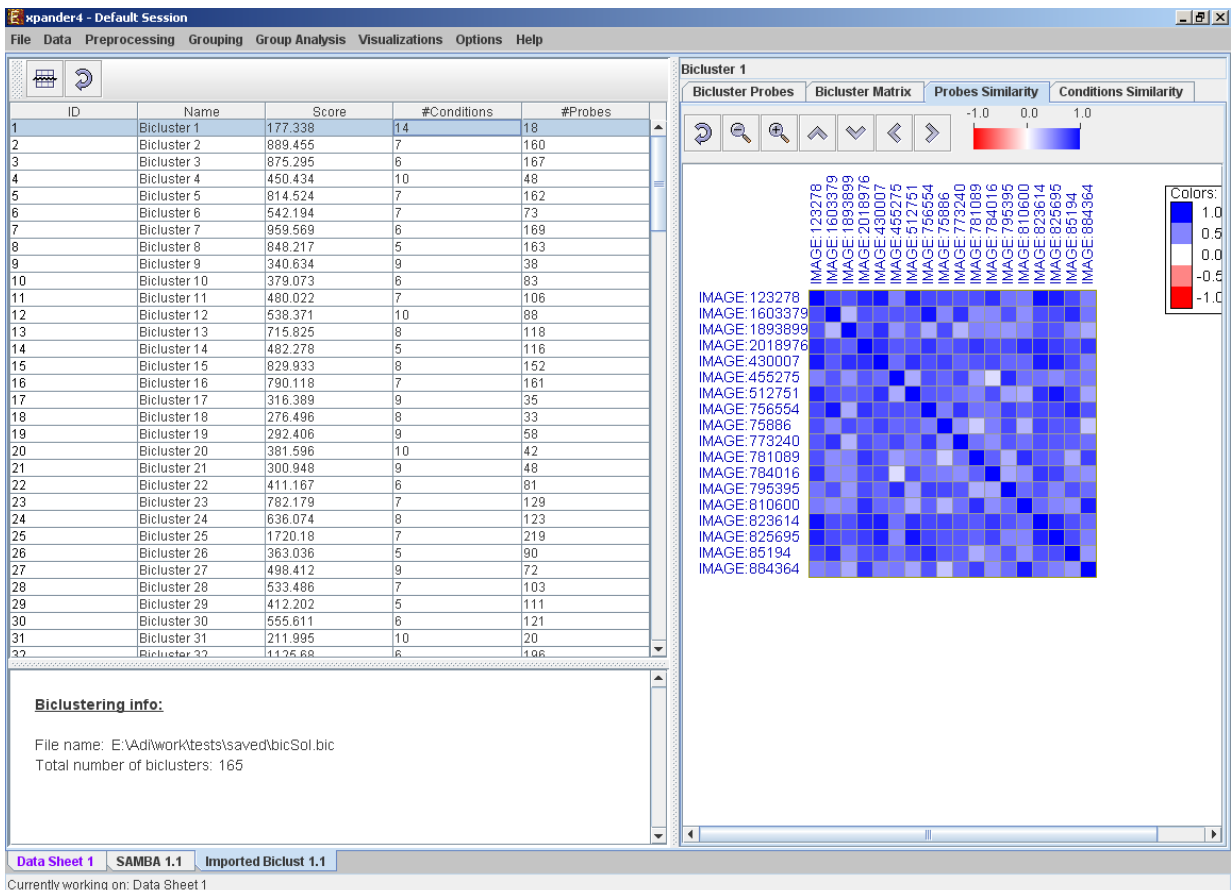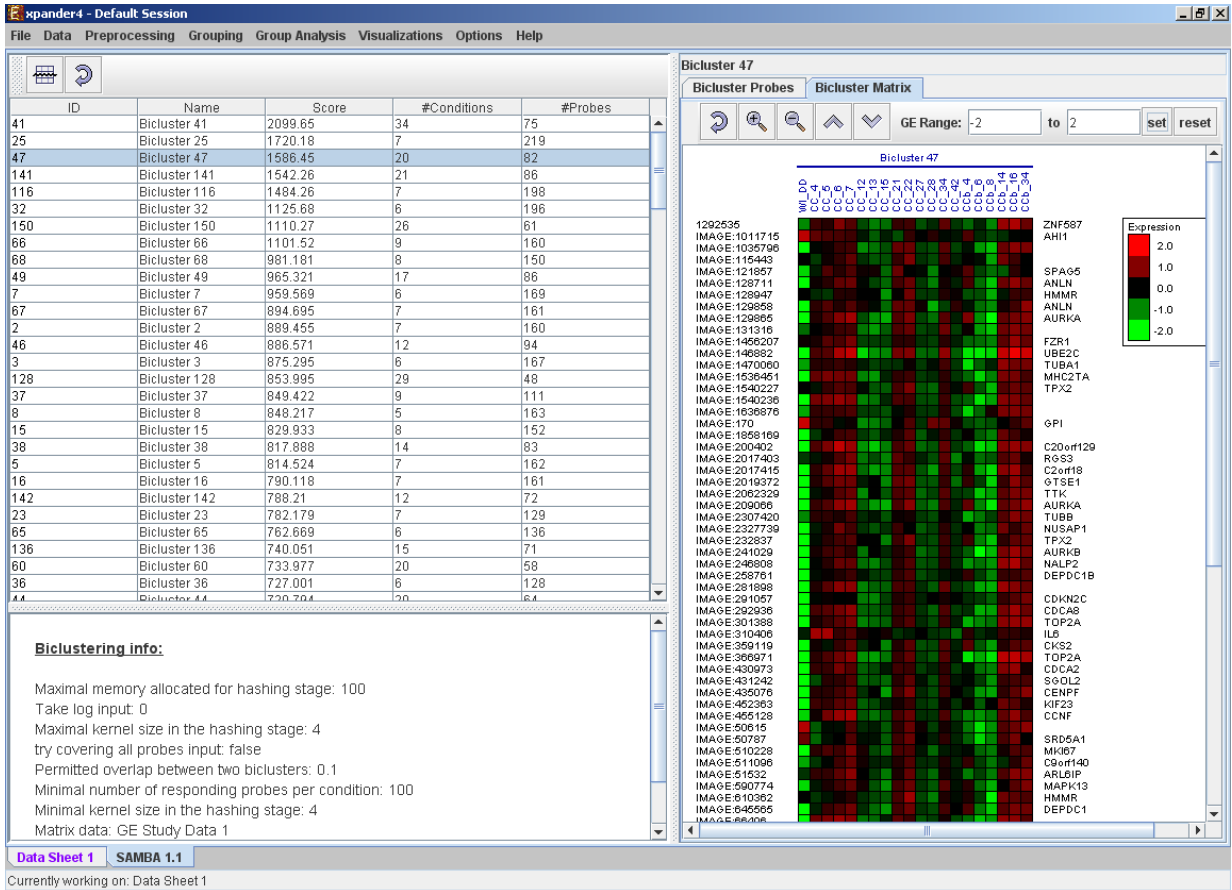

After performing group analysis (for details see the [Group Analysis Tools](#) section), if enrichment has been detected in the selected bicluster, the corresponding histogram and analysis information are added to the single bicluster view, and a column is added to the expression matrix display for each enrichment class, stating for each probe, whether it belongs to that class.

A biclustering solution can be saved using the *Grouping >> Bi-Clustering >> Save Solution*, and reloaded using the *Grouping >> Bi-Clustering >> Load Solution*. For a format of the solution file, please refer to the [File Formats](#) section:

Network Based Grouping of GE Data


The goal here is to detect groups of genes that demonstrate similar expression patterns and are also highly connected in a given interactions network.

Before operating this tool, an interactions network in .SIF format must be loaded. This can be done by selecting *Data>>Load network*.

Network based grouping is performed by Expander using the Matisse algorithm (for details see the [References](#) section). The groups detected by Matisse are referred to as “modules” and may contain also genes that exist in the network, but are not present in the filtered GE data (referred to as “Back nodes”).

To use the more advanced, stand-alone version of MATISSE (with higher flexibility), please refer to [the Matisse home page](#).

In order to apply the Matisse algorithm to the data select *Grouping>>Network >>Matisse*. The following dialog box will appear:


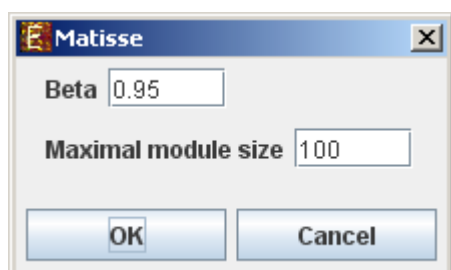


It enables the configuration of some of the parameters for the algorithm:

Field	Description
Beta	The fraction of gene pairs that are expected to be strongly co-expressed in each module
Maximal module size	The maximum size for a detected module.

Upon clicking ‘OK’ in the dialog box, the Matisse algorithm is operated on the dataset.

After running the algorithm, the solution is displayed in a new tab, which is added to the main window. The view is similar to the clustering results display. In the display, back nodes (genes that appear in the network, but not in the GE data) are marked in yellow.

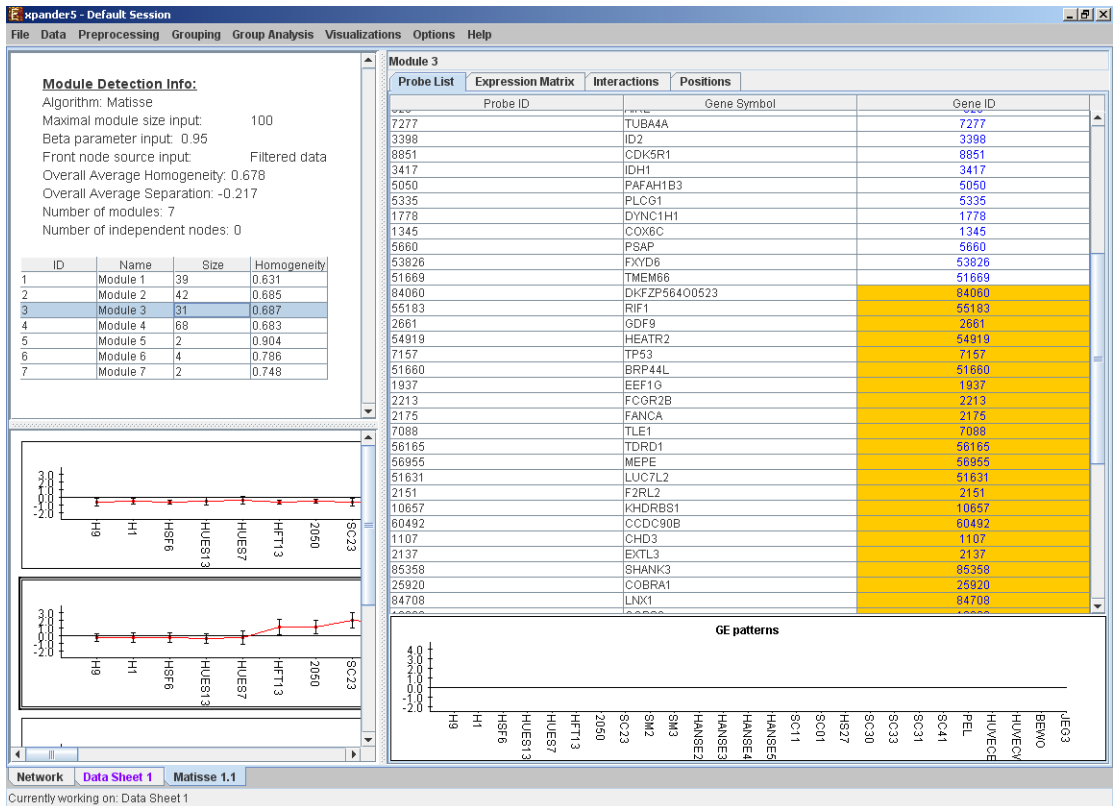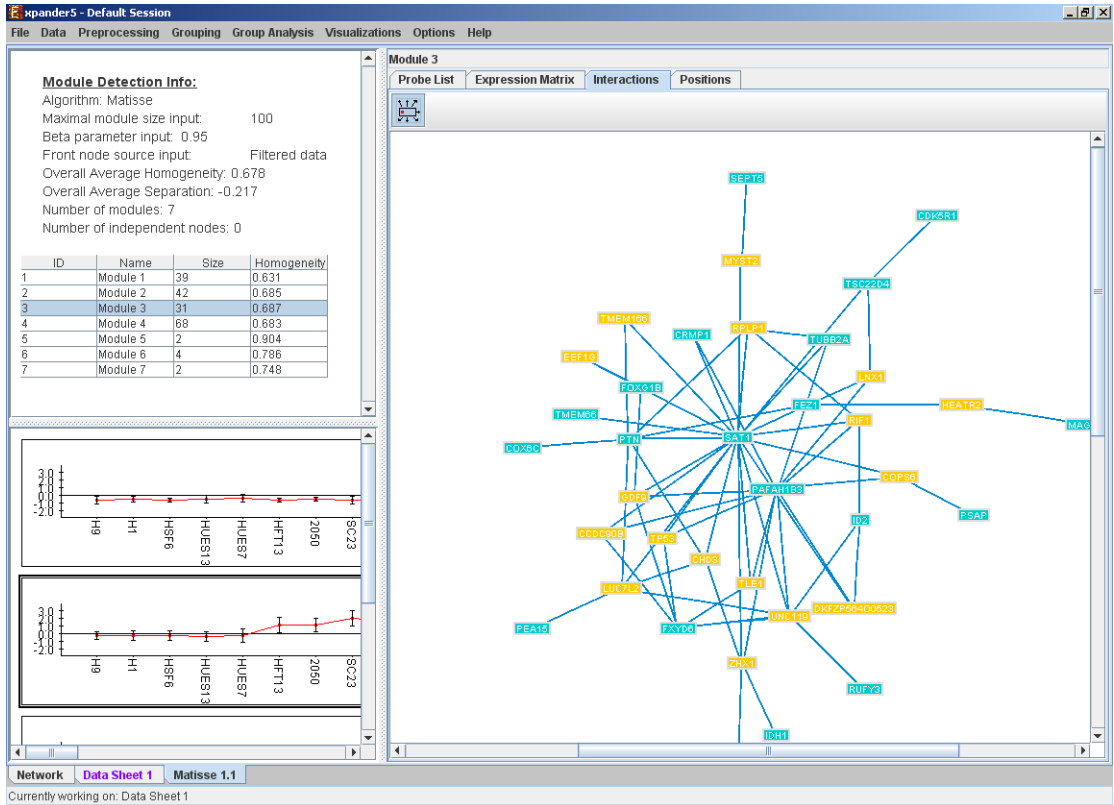

After performing group analysis (for details see the [Group Analysis Tools](#) section), if enrichment has been detected in the selected module, the corresponding histogram and analysis information are added to the single module view, and a column is added to the expression matrix display for each enrichment class, stating for each probe, whether it belongs to that class.

A network-based grouping solution can be saved using the *Grouping >> Network >> Save Solution*, and reloaded using the *Grouping >> Network >> Load Solution*. For a format of the solution file, please refer to the [File Formats](#) section:

Group Analysis Tools

The following analysis can be performed on gene sets, clusters, biclusters, network based modules, similarity based groups, or the filtered dataset (the analyzed set of probes as one set). Before operating any of the group analysis operation (not including the "General enrichment analysis"), the data files for the relevant organism should be downloaded. The files for a specific organism are supplied as one single zip file that you need to download from the Expander download page, section "Organism specific data" (can be reached from inside Expander: Help >> Open download page, or from updates.html file in your Expander directory). Then, you extract the file into "Expander/organisms" directory. The relevant directories will be built by the extraction, inside your "Expander/organisms" directory. For example, after extracting human data, you should have "Expander/organisms/human" directory.

Functional Analysis

This tool performs basic statistical analysis on the distribution of functions of genes within each cluster. The functions of the genes are determined according to annotation files (GO), which can be downloaded from the EXPANDER download page (see the [Supplied Files](#) section). To perform this analysis, Expander utilizes the TANGO software, which performs hyper-geometric enrichment tests and corrects for multiple testing by bootstrapping and estimating the empirical p-value distribution for the evaluated sets.

Before operating functional analysis the annotation files for the relevant organism should be downloaded from the download page (more details at introduction of [Group Analysis Tools](#)). To perform the analysis, select *Group Analysis >> Functional Analysis >> TANGO*. The following dialog box will appear:

The following table specifies the different parameters that can be set via this dialog box:


Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Include back nodes	Include genes that are part of the module' but not included in the GE data (Relevant only if the analysis is performed on modules, detected by network based algorithm)
Focus on	Can be used to select annotation subtypes that are of interest (Process, Function and Location). And the analysis will focus on these types only.
Ignore classes over the size of	This parameter states the level in the GO tree at which annotations are too general (class size indicates how general it is) and are thus no longer interesting.
Number of iterations in algorithm	The number of random sampling performed by the algorithm. Increasing this parameter, will increase runtime and will provide higher resolution on corrected p-Values. I.e., corrected p-Values will range between $1/\text{#iterations}$ and 1.
Background set	Determines the set of genes that will be used as background in the analysis. Options are: all genes (of the relevant organism), original input data, filtered data or background set from file (see the

	Files Format section for details regarding the format of an external background set).
Corrected p-value threshold	A functional class will be considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than this threshold. The value in this field should be at least 1/1000, since the TANGO algorithm performs 1000 bootstraps in order to estimate the corrected p-value.

Upon clicking 'OK' in the dialog box, the TANGO algorithm is operated.

After functional analysis is performed a functional analysis solution visualization tab is added to the main window. It contains general information regarding the analysis, a sort-able table holding all detected enrichments (set ID, functional class, p-value, etc.) and a multi-histogram panel along with a color index (mapping each color to a corresponding functional class). The multi-histogram panel contains one histogram for each probe/gene set/group in which enrichment has been detected. Each histogram contains a column for each significant (more frequent than would be expected by random) functional class. The definition of significant depends on the user's selection of threshold p-value i.e., a functional class is considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than the preset threshold p-value.

The height of the column is proportional to the significance of this enrichment (i.e. height = $-\log(\text{raw p-value})$). The frequency in set (frequency of genes of a functional class within the examined set, in %) is written on top of the column. Upon clicking on a column, a dialog box is displayed containing the class name, raw p-value, corrected p-value, and a list of the genes in the cluster/bi-cluster that belong to the class. Upon clicking on one of the gene Ids in the table, a relevant web page with information regarding this gene is displayed. The display tool tip shows the cluster number, size and homogeneity.

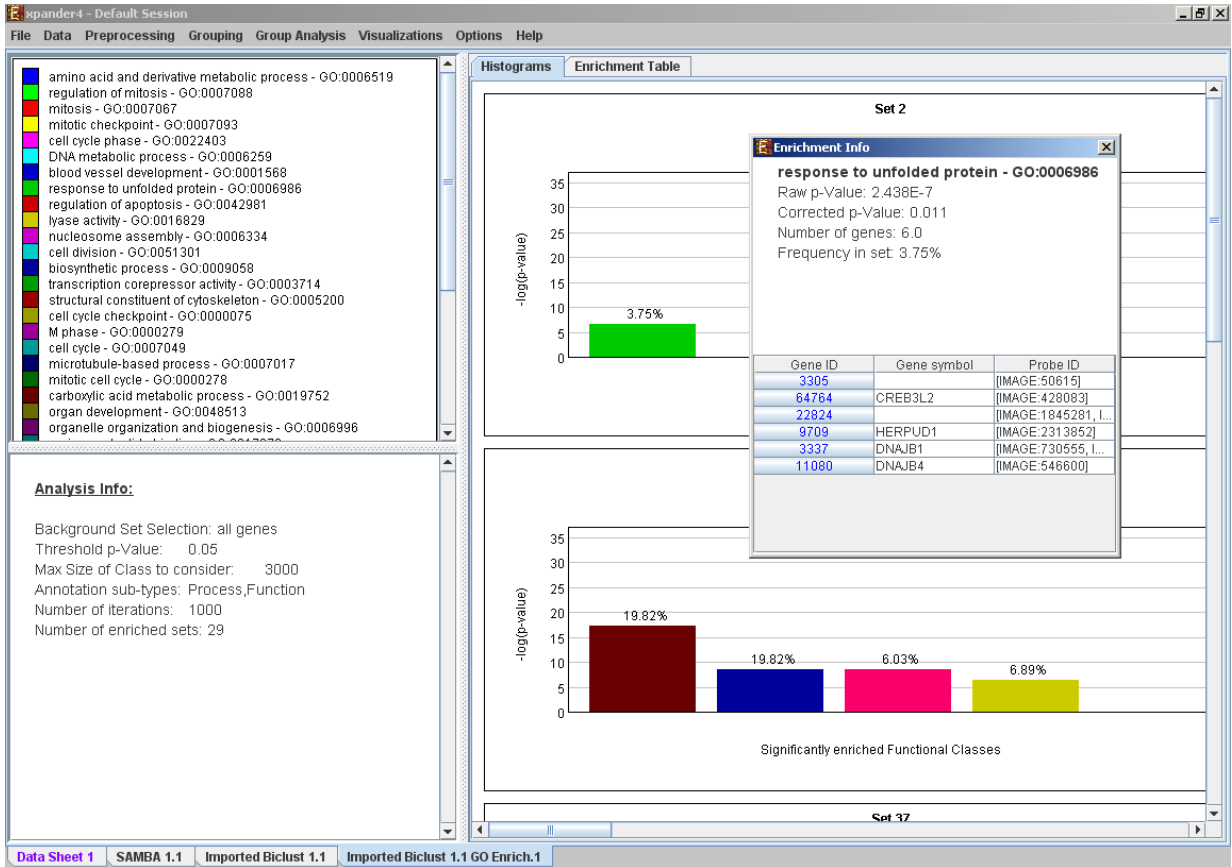

Expander4 - Default Session

File Data Preprocessing Grouping Group Analysis Visualizations Options Help

- amino acid and derivative metabolic process - GO:0006519
- regulation of mitosis - GO:0007088
- mitosis - GO:0007067
- mitotic checkpoint - GO:0007093
- cell cycle phase - GO:0022403
- DNA metabolic process - GO:0006259
- blood vessel development - GO:0001568
- response to unfolded protein - GO:0006986
- regulation of apoptosis - GO:0042981
- lyase activity - GO:0016829
- nucleosome assembly - GO:0006334
- cell division - GO:0051301
- biosynthetic process - GO:0009058
- transcription corepressor activity - GO:0003714
- structural constituent of cytoskeleton - GO:0005200
- cell cycle checkpoint - GO:0000075
- M phase - GO:0000279
- cell cycle - GO:0007049
- microtubule-based process - GO:0007017
- mitotic cell cycle - GO:0000278
- carboxylic acid metabolic process - GO:0019752
- organ development - GO:0048513
- organelle organization and biogenesis - GO:0006996

Analysis Info:

Background Set Selection: all genes
 Threshold p-Value: 0.05
 Max Size of Class to consider: 3000
 Annotation sub-types: Process,Function
 Number of iterations: 1000
 Number of enriched sets: 29

Histograms Enrichment Table

Set	Enriched with	#genes	p-value	Frequency in set (%)
2	response to unfolded protein - GO:0...	6	2.438E-7	3.75
6	response to unfolded protein - GO:0...	5	1.97E-7	6.84
7	response to unfolded protein - GO:0...	8	3.045E-10	4.73
8	response to unfolded protein - GO:0...	7	8.562E-9	4.29
13	amino acid and derivative metabolic...	10	6.185E-7	8.47
13	response to unfolded protein - GO:0...	5	1.225E-6	4.23
14	carboxylic acid metabolic process - ...	23	4.637E-18	19.82
14	biosynthetic process - GO:0009058	23	2.47E-9	19.82
14	ligase activity, forming carbon-oxyg...	7	2.604E-9	6.03
14	lyase activity - GO:0016829	8	3.57E-7	6.89
15	regulation of apoptosis - GO:0042981	15	6.065E-8	9.86
23	transcription corepressor activity - G...	7	4.006E-7	5.42
24	cell cycle arrest - GO:0007050	6	8.751E-7	4.87
36	organ development - GO:0048513	18	1.629E-7	14.06
36	organ morphogenesis - GO:0009887	10	7.85E-7	7.81
37	mitotic cell cycle - GO:0000278	21	1.222E-22	18.91
37	mitosis - GO:0007067	19	1.783E-22	17.11
37	cell cycle - GO:0007049	28	3.417E-20	25.22
37	cell division - GO:0051301	17	4.073E-18	15.31
37	microtubule-based process - GO:00...	12	3.442E-12	10.81
37	spindle organization and biogenesi...	6	1.539E-10	5.4
37	cell cycle checkpoint - GO:0000075	5	1.202E-6	4.5
38	organ development - GO:0048513	20	1.298E-11	24.09
38	response to wounding - GO:0009611	12	9.471E-10	14.45
38	wound healing - GO:0042060	7	2.943E-8	8.43
38	cell motility - GO:0006928	10	6.856E-8	12.04
38	anatomical structure formation - GO:...	7	2.706E-7	8.43
38	blood vessel development - GO:000...	7	2.706E-7	8.43
39	structural constituent of cytoskeleton...	6	1.828E-7	7.31
40	cell motility - GO:0006928	9	1.17E-7	15.0
41	M phase - GO:0000279	27	1.241E-37	36.0
41	mitosis - GO:0007067	25	1.567E-36	33.33
41	mitotic cell cycle - GO:0000278	27	1.942E-36	36.0
41	cell cycle - GO:0007049	32	3.511E-30	42.66
41	cell division - GO:0051301	22	4.338E-29	29.33
41	microtubule-based process - GO:00...	12	9.605E-14	16.0
41	spindle organization and biogenesi...	7	1.548E-13	9.33
41	cell cycle checkpoint - GO:0000075	6	5.662E-9	8.0
42	nucleosome assembly - GO:0006334	6	1.76E-8	10.52
42	DNA metabolic process - GO:00062...	11	2.027E-7	19.29
44	nucleosome assembly - GO:0006334	7	3.393E-10	10.93
44	DNA metabolic process - GO:00062...	12	2.0E-8	18.75
44	cellular component organization an...	19	2.214E-7	29.68
44	organelle organization and biogene...	12	1.318E-6	18.75

Data Sheet 1 SAMBA 1.1 Imported Biclust 1.1 Imported Biclust 1.1 GO Enrich.1

Currently working on: Data Sheet 1

Annotation files are currently supplied with EXPANDER for yeast, human, mouse, rat, fly, zebrafish, c-elegans, Arabidopsis, chicken and E. coli, and are updated on a regular basis (for more information, refer to the [Supplied Files](#) section).

The results of this analysis can be exported to a text file (using the display menu item *Group Analysis >> Functional Analysis >> Export Results*).

Promoter Analysis

This tool identifies TFs whose binding sites are significantly over-represented in a given set of promoters (i.e. cluster or bicluster). To perform this analysis Expander utilizes the PRIMA (PRomoter Integration in Microarray Analysis) software which performs a statistical analysis on the distribution of transcription factor motifs in the promoters of genes within each cluster or bicluster. To achieve this, PRIMA uses preprocessed TF fingerprint files, which can be downloaded from the EXPANDER download-page (see the [Supplied Files](#) section), and are updated on a regular basis. For details regarding the PRIMA software see the [References](#) section.

Before operating promoter analysis, the TF fingerprint file for the relevant organism should be downloaded from the download page (more details at introduction of [Group Analysis Tools](#)). To perform the analysis select *Group Analysis >> Promoter Analysis >> PRIMA*. The following dialog box will appear:

Promoter Analysis

Perform Analysis On: CLICK 1.1

Fingerprints file: fp.Masked.fp1000.txt

PWM file: expander_pwmsIndSet_0.25_Hs_v27.l1200_mammals_298.txt

Promoter sequences file: Hs_Proms_v27.all.txt

Consider hits from -1000 to 200 Ignore coding regions

Background set

All Genes

Original Data

Filtered Data

From File Browse

p-value threshold: 1.0E-4 Multiple tests correction: None

Save results as: Browse

OK Cancel

The following table specifies the different parameters that can be set via this dialog box:


Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Include back nodes	Include genes that are part of the module' but not included in the GE data (Relevant only if the analysis is performed on modules, detected by network based algorithm)
Fingerprints file	Automatically set according to the selection of the organism.
PWM file	Automatically set according to the selection of the organism.
Promoter sequences file	Contains the gene sequences that are used for the TF binding sites display. Automatically set according to the selection of the organism.
Hits range	Determines which regions of the gene are to be analyzed. The possible range depends on the investigated organism (i.e. on the information provided in the TF fingerprint files), and is specified in the Supplied Files section.
Background set	Determines the set of genes that will be used as background in the analysis. Options are: all genes (of the relevant organism), original input data, filtered data or background set from file (see the Files Format section for details regarding the format of an external background set).
Threshold p-value	A TF's binding site will be considered significantly enriched in a cluster if its corrected p-value is lower than this threshold.
Multiple tests correction	Can be set to Bonferroni or None (when set to Bonferroni the corrected p-values are the ones that are compared to the threshold p-value).
Save results as	When filled, the program results are saved in stated txt file.

After promoter analysis is performed, a promoter analysis solution visualization tab is added to the main window. It contains general information regarding the analysis, a sort-able table holding all detected enrichments (set ID, TF binding site, p-value, etc.) and a multi-histogram panel along with a color index (mapping each color to a corresponding TF binding site). The multi-histogram panel contains one histogram for each probe/gene set/group in which enrichment has been detected. Each histogram contains a column for each significant (more frequent than would be expected by random) TF binding site. The definition of significant depends on the user's selection of threshold p-value. i.e., a TF binding site is considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than the preset threshold p-value.

The height of a column is proportional to the significance of this enrichment (i.e. height = $-\log(p\text{-value})$), and the frequency ratio (frequency in set divided by frequency

in background) is written on top of the column. Upon clicking on a column, a dialog box is displayed containing:

TF accession number in TRANSFAC DB [TF name], p-value, % of covered promoters in cluster, relative frequency (frequency in cluster divided by frequency in background set) and a list of the genes in the cluster which contain the motif in their promoters. Upon clicking on one of the gene Ids in the table, a relevant web page with information regarding this gene is displayed. The display tool tip shows the cluster number, size and homogeneity.

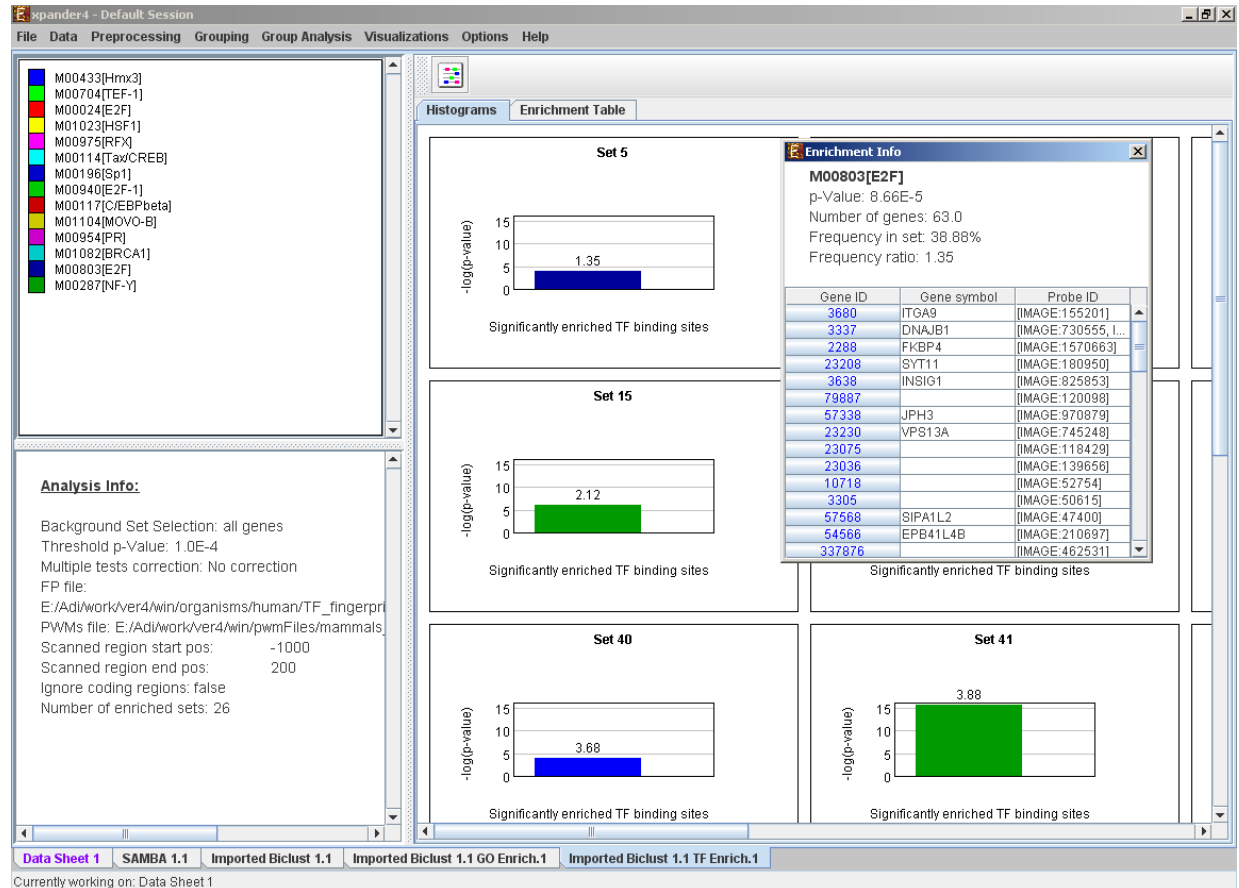

Expander4 - Default Session

File Data Preprocessing Grouping Group Analysis Visualizations Options Help

Color Index:

- M00433[Hmx3]
- M00704[TEF-1]
- M00024[E2F]
- M01023[HSF1]
- M00975[RFx]
- M00114[Taw/CREB]
- M00196[Sp1]
- M00940[E2F-1]
- M00117[C/EBPbeta]
- M01104[MOVO-B]
- M00954[PR]
- M01082[BRCA1]
- M00803[E2F]
- M00287[NF-Y]

Analysis Info:

Background Set Selection: all genes
Threshold p-Value: 1.0E-4
Multiple tests correction: No correction
FP file:
E:/Adi/work/ver4/win/organisms/hum
PWMs file: E:/Adi/work/ver4/win/pwm
Scanned region start pos: -1
Scanned region end pos: 20
Ignore coding regions: false
Number of enriched sets: 26

Enrichment Table

Set	Enriched with	#genes	p-value	Enrichment factor
5	M00803[E2F]	63	8.66E-5	1.352
6	M01023[HSF1]	14	2.25E-5	3.001
7	M01023[HSF1]	22	8.07E-6	2.201
7	M00196[Sp1]	76	9.6E-8	1.418
7	M01104[MOVO-B]	55	4.94E-5	1.533
8	M01023[HSF1]	22	1.21E-5	2.277
8	M00287[NF-Y]	34	2.25E-5	1.723
13	M00117[C/EBPbeta]	15	5.22E-5	2.952
13	M00287[NF-Y]	31	3.52E-8	2.277
15	M00287[NF-Y]	39	7.38E-7	2.122
19	M00287[NF-Y]	17	3.22E-5	2.775
22	M00287[NF-Y]	21	6.7E-5	2.286
23	M00287[NF-Y]	27	2.1E-5	1.984
37	M00287[NF-Y]	29	8.28E-7	2.435
40	M00433[Hmx3]	14	6.76E-5	3.682
41	M00287[NF-Y]	35	2.49E-16	3.882
42	M00024[E2F]	9	3.78E-5	4.147
42	M00940[E2F-1]	11	6.95E-5	3.208
42	M00287[NF-Y]	22	2.28E-7	2.939
44	M00287[NF-Y]	23	9.05E-8	3.072
46	M00287[NF-Y]	34	4.74E-11	3.122
47	M00287[NF-Y]	36	1.78E-15	3.527
52	M01082[BRCA1]	8	4.11E-5	4.219
52	M00287[NF-Y]	16	4.83E-8	3.359
53	M00803[E2F]	39	4.49E-5	1.546
79	M00954[PR]	13	6.32E-5	3.926
90	M00975[RFx]	7	7.2E-6	4.668
96	M00704[TEF-1]	15	5.91E-6	2.349
128	M00287[NF-Y]	15	6.35E-6	2.593
131	M00287[NF-Y]	16	4.47E-9	3.617
134	M00287[NF-Y]	20	9.05E-10	4.054
136	M00114[Taw/CREB]	16	1.09E-5	2.571
141	M00287[NF-Y]	36	1.04E-15	3.469


Currently working on: Data Sheet 1

יזגה בגן הרושטלר - Mozilla Firefox

After performing promoter analysis, TF binding sites can be viewed by selecting *Group Analysis >> Promoter Analysis >> View Binding Sites* OR by pressing the



After selecting the gene group (cluster/bi-cluster etc.) to be viewed, a separate frame is displayed, containing a line to represent each of the genes in the group, and a colored rectangle, to represent each binding site. A color index appears on the right, mapping each color to the corresponding TF (PWM). A check box next to each of the entries in the color index allows hiding any of the PWMs, and a radio button next to each of the entries in the color index allows sorting the genes in the display according to the number of hits of the corresponding TF. The toolbar contains tools for vertical and horizontal zooming. If a sequence file had been selected via the promoter analysis input dialog, the actual sequence will be displayed when the zoom factor (scale) allows it.

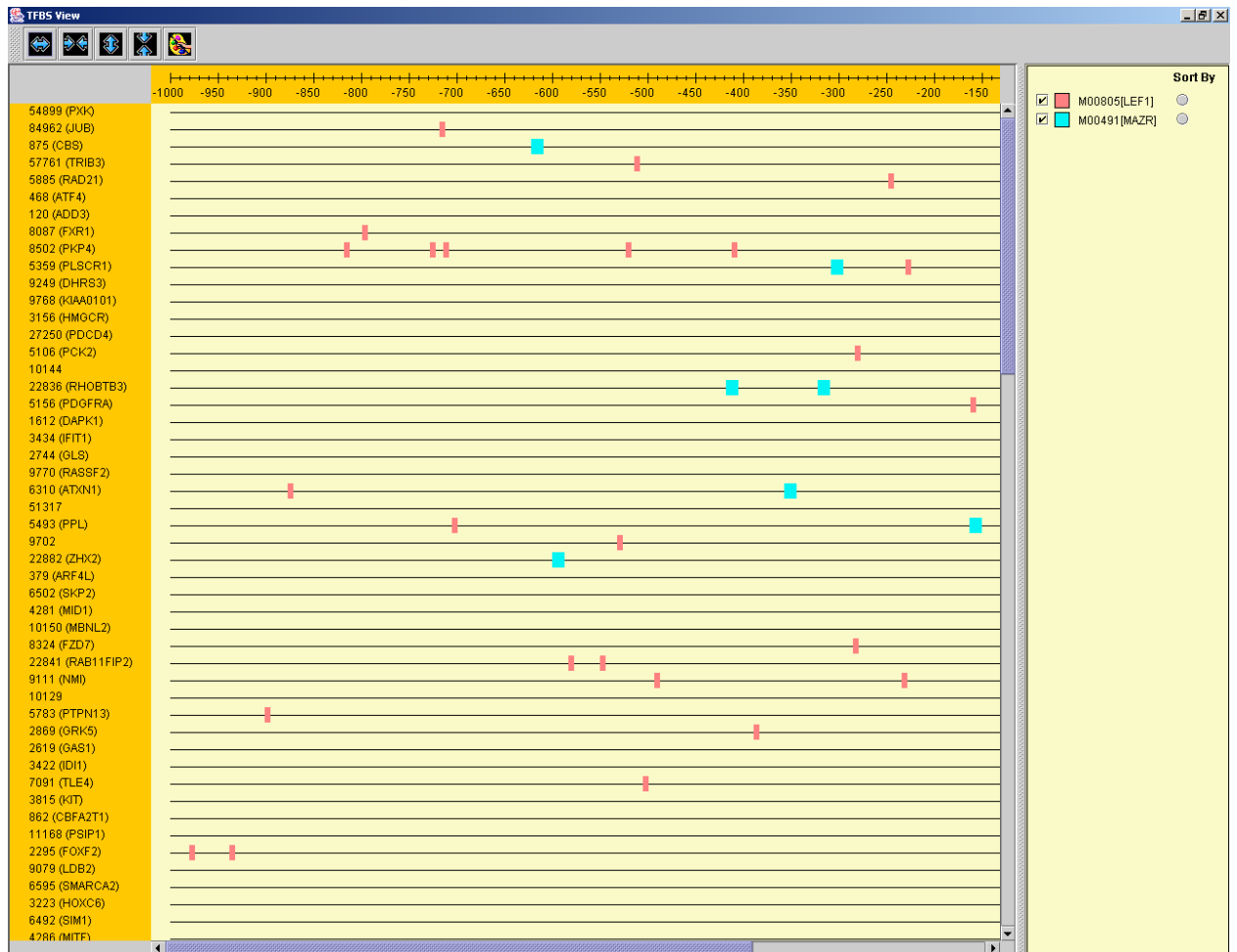

TF motif fingerprint files and promoter sequence files are currently supplied with EXPANDER for yeast, human, mouse, rat, fly, zebrafish, c-elegans, arabidopsis and chicken, and are updated on a regular basis (for more information, refer to the [Supplied Files](#) section).


The results of this analysis can be exported to a text file (using the display menu item *Group Analysis >> Promoter Analysis >> Export Results*).

Location Enrichment Analysis

This tool performs basic statistical analysis on the distribution of chromosomal locations of genes within each group. The locations of the genes are specified in organism-specific data files, which can be downloaded from the EXPANDER download-page (see the [Supplied Files](#) section).

Before operating location analysis, the location data for the relevant organism should be downloaded from the download page (more details at introduction of [Group Analysis Tools](#)). In this analysis, hyper-geometric enrichment tests are performed, and the results can be (if requested) corrected for multiple testing using the Bonferroni correction.

To perform the analysis, select *Group Analysis >> Location Analysis >> Detect Enrichment*. The following dialog box will appear:

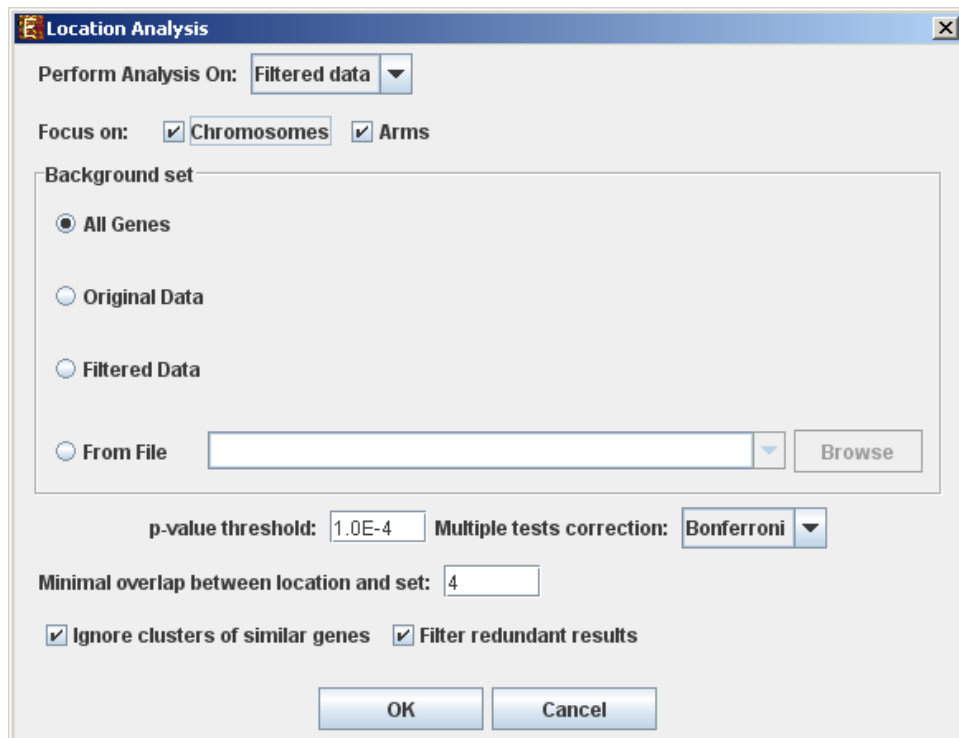

The following table specifies the different parameters that can be set via this dialog box:


Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Include back nodes	Include genes that are part of the module' but not included in the GE data (Relevant only if the analysis is performed on modules, detected by network based algorithm)
Focus on (Chromosomes, Arms*, Bands*)	Location types to perform analysis on.
Background set	Determines the set of genes that will be used as background in the analysis. Options are: all genes (of the relevant organism), original input data, filtered data or background set from file (see the Files Format section for details regarding the format of an external background set).
p-value threshold	A category/attribute will be considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than this threshold.
Multiple tests correction	Can be set to Bonferroni or None (when set to Bonferroni the corrected p-values are the ones that are compared to the threshold p-value).
Minimal overlap between	The minimal number of genes from a group (cluster/

category and set	bi-cluster/module etc.) expected to be categorized/attributed by an attribute in order for its enrichment to be accepted.
Ignore clusters of similar genes*	If selected, genes from known homology clusters are not included in the analysis.
Filter redundant results	If selected, the results are filtered, so that out of two enrichments of overlapping areas in the same group, only one is selected (the most significant one).

* If relevant data exists

After the analysis is performed an enrichment analysis solution visualization tab is added to the main window. It contains general information regarding the analysis, a sort-able table holding all detected enrichments (set ID, enrichment category, p-value, etc.) and a multi-histogram panel along with a color index (mapping each color to a corresponding location). The multi-histogram panel contains one histogram for each probe/gene group in which enrichment has been detected. Each histogram contains a column for each significant (more frequent than would be expected by random) location. The definition of significant depends on the user's selection of threshold p-value i.e., a category is considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than the preset threshold p-value.


The height of the column is proportional to the significance of this enrichment (i.e. height = $-\log(\text{raw p-value})$), and the frequency ratio (frequency in set divided by frequency in background) is written on top of the column. Upon clicking on a column, a dialog box is displayed containing the location, corrected p-value, and a list of the genes in the group that are mapped to this location. Upon clicking on one of the gene IDs in the table, a relevant web page with information regarding this gene is displayed.

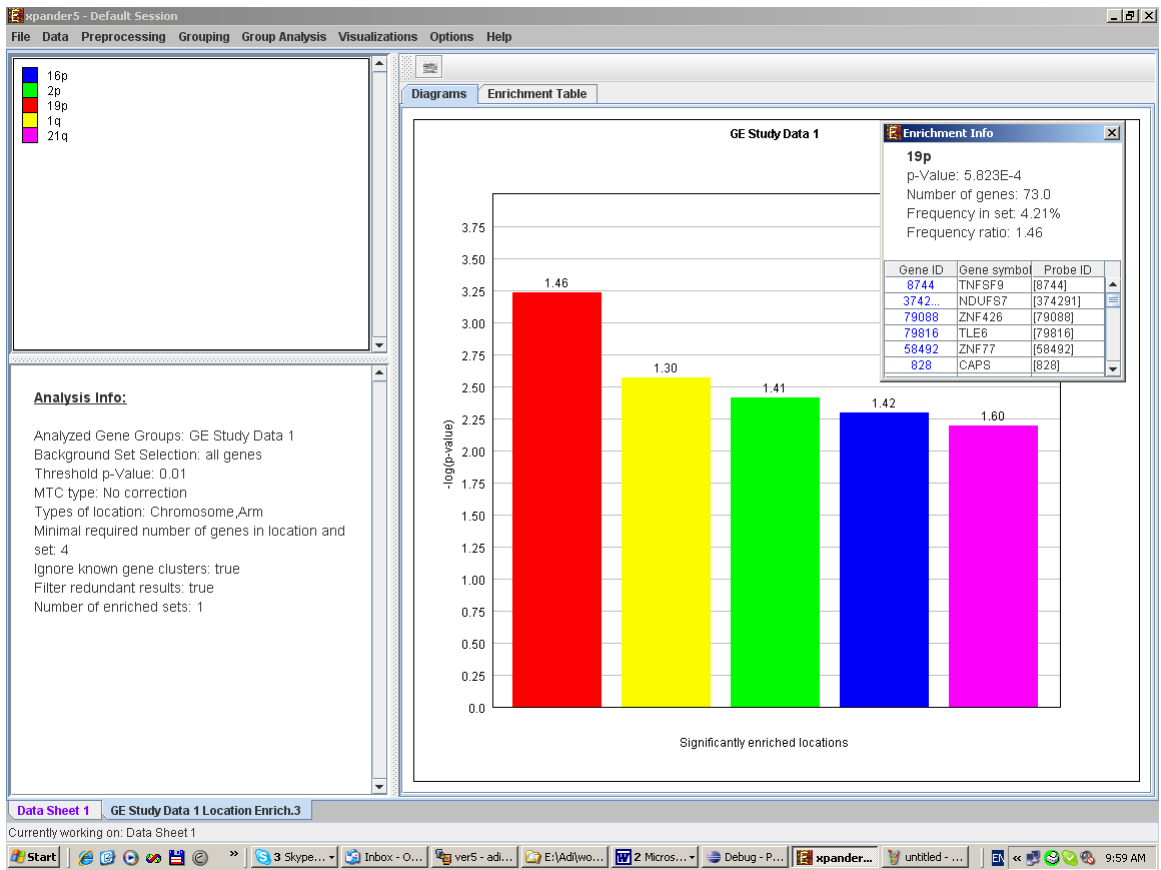


Set	Enriched with	#genes	p-Value	Enrichment factor
GE Study Data 1	1q	106	0.0020	1.302
GE Study Data 1	2p	59	0.0030	1.417
GE Study Data 1	16p	53	0.0050	1.429
GE Study Data 1	21q	30	0.0060	1.606
GE Study Data 1	19p	73	5.823E-4	1.464

Analysis Info:

Analyzed Gene Groups: GE Study Data 1
 Background Set Selection: all genes
 Threshold p-Value: 0.01
 MTC type: No correction
 Types of location: Chromosome_Arm
 Minimal required number of genes in location and set: 4
 Ignore known gene clusters: true
 Filter redundant results: true
 Number of enriched sets: 1

After performing location enrichment analysis, the locations can be viewed by selecting *Group Analysis* >> *Location Analysis* >> *View Locations* OR by pressing the toolbar button (). After selecting the gene group (cluster/bi-cluster etc.) to be

viewed, a separate frame is displayed, containing an image of all chromosomes on which the positions of the genes in the group are marked. If the gene is located on an area that was identified to be enriched in that group, its position is marked in the same color to this area the enrichment results histogram.

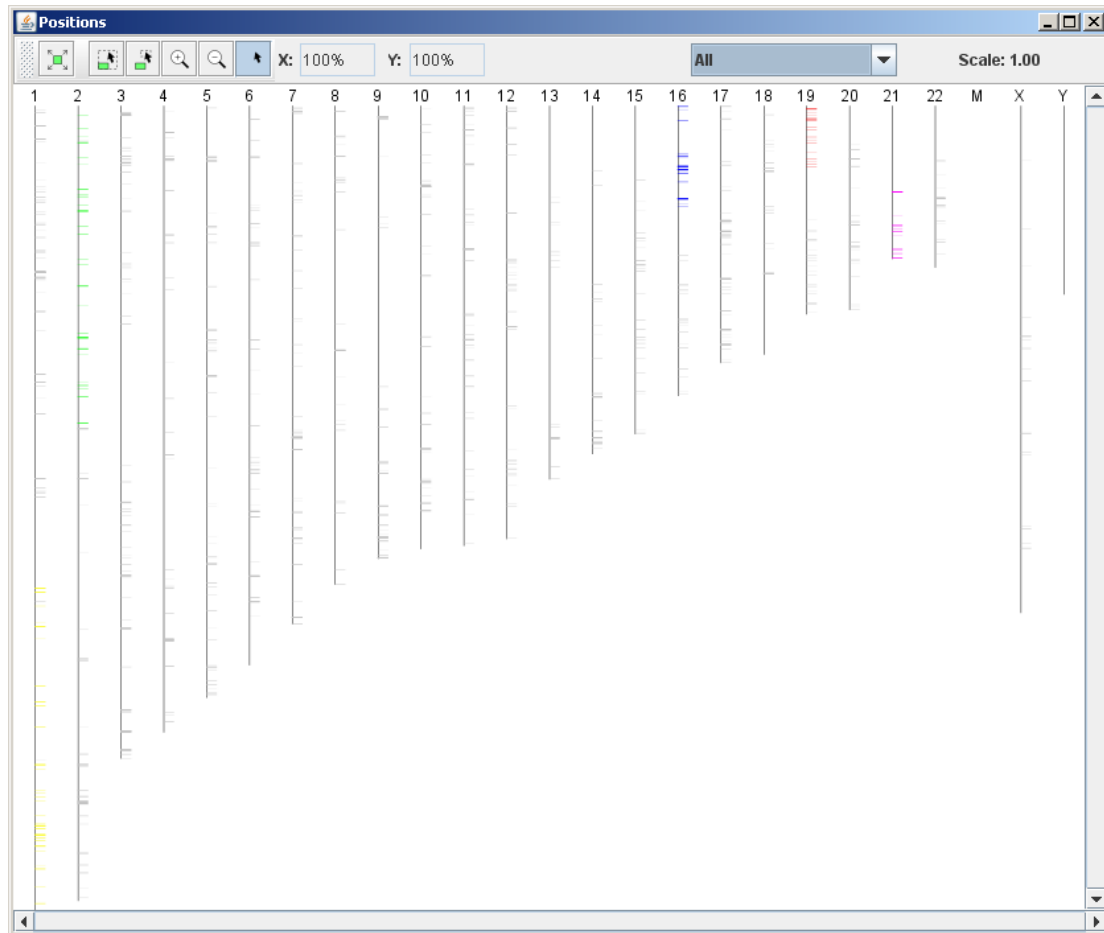

The results of this analysis can be exported to a text file (using the display menu item *Group Analysis >> Location Analysis >> Export Results*).

miRNA Targets Enrichment Analysis

This tool performs a statistical analysis on the distribution of miRNA target gene within each group. The miRNA targets information is supplied in organism-specific data files, which can be downloaded from the EXPANDER download-page (see the [Supplied Files](#) section). For this analysis, Expander utilizes the FAME algorithm, which performs empirical tests using a sampling technique (random permutations) to estimate the empirical p-value distribution for the evaluated groups. This is done while accounting for biases in the 3' UTR sequences

Before operating miRNA enrichment analysis, the location data for the relevant organism should be downloaded from the download page (more details at introduction of [Group Analysis Tools](#)). In this analysis, hyper-geometric enrichment tests are performed, and the results can be (if requested) corrected for multiple testing using the Bonferroni correction.

To perform the analysis, select *Group Analysis >> miRNA Analysis >> FAME*. The following dialog box will appear:


The following table specifies the different parameters that can be set via this dialog box:

Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Include back nodes	Include genes that are part of the module' but not included in the GE data (Relevant only if the analysis is performed on modules, detected by network based algorithm)
Enrichment Direction	Allows to choose between searching for over-represented targets and searching for under-represented targets.
Use context scores	If context scores are used, FAME will assign a higher weight to miRNA-gene pairs for which at least one target site has a high maximal context score (see References section for further details).
Number of Iterations	The number of random permutations used for the empirical tests.
Background set	Determines the set of genes that will be used as background in the analysis. Options are: all genes (of the relevant organism), original input data,

	filtered data or background set from file (see the Files Format section for details regarding the format of an external background set).
p-value threshold	A category/attribute will be considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than this threshold.
Multiple tests correction	Can be set to FDR, Bonferroni or None (when set to Bonferroni/FDR the corrected p-values are the ones that are compared to the threshold p-value).
Minimal overlap between targets and group	The minimal number of genes from a group (cluster/ bi-cluster/module etc.) expected to be categorized/attributed by an attribute in order for its enrichment to be accepted.

After the analysis is performed an enrichment analysis solution visualization tab is added to the main window. It contains general information regarding the analysis, a sort-able table holding all detected enrichments (group name, enriched miRNA target, p-value, etc.) and a multi-histogram panel along with a color index (mapping each color to a corresponding miRNA). The multi-histogram panel contains one histogram for each probe/gene group in which enrichment has been detected. Each histogram contains a column for each significant (more frequent than would be expected by random) miRNA target. The definition of significant depends on the user's selection of threshold p-value i.e., an mRNA target is considered significantly enriched in a group of genes if its corrected p-value is lower than the selected threshold p-value.

The height of the column is proportional to the significance of this enrichment (i.e. height = $-\log(\text{raw p-value})$), and the frequency ratio (frequency in set divided by frequency in background) is written on top of the column. Upon clicking on a column, a dialog box is displayed containing the miRNA name, corrected p-value, and a list of the genes in the group that are mapped to this location. Upon clicking on one of the gene ids in the table, a relevant web page with information regarding this gene is displayed.

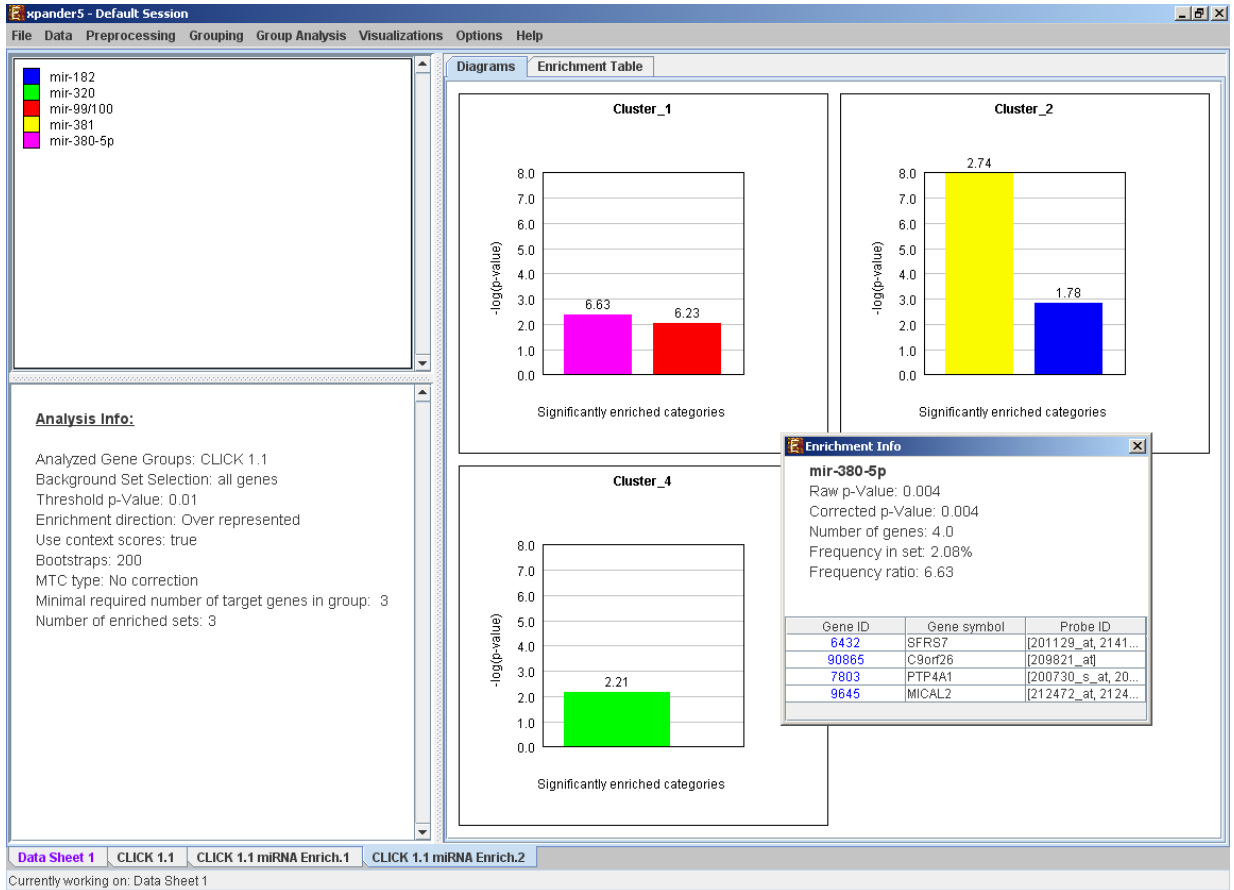

xpander5 - Default Session

File Data Preprocessing Grouping Group Analysis Visualizations Options Help

Legend:

- mir-182
- mir-320
- mir-99/100
- mir-381
- mir-380-5p

Diagrams Enrichment Table

Set	Enriched with	#genes	p-Value	Enrichment factor
Cluster_1	mir-380-5p	4	0.0040	6.637
Cluster_1	mir-99/100	3	0.0090	6.239
Cluster_2	mir-381	26	1.045E-8	2.748
Cluster_2	mir-182	17	0.0010	1.788
Cluster_4	mir-320	8	0.0060	2.211

Analysis Info:

Analyzed Gene Groups: CLICK 1.1
 Background Set Selection: all genes
 Threshold p-Value: 0.01
 Enrichment direction: Over represented
 Use context scores: true
 Bootstraps: 200
 MTC type: No correction
 Minimal required number of target genes in group: 3
 Number of enriched sets: 3

Data Sheet 1 CLICK 1.1 CLICK 1.1 miRNA Enrich.1 CLICK 1.1 miRNA Enrich.2

Currently working on: Data Sheet 1

The results of this analysis can be exported to a text file (using the display menu item *Group Analysis >> miRNA Analysis >> Export Results*).

General Enrichment Analysis

This tool performs basic statistical analysis on the distribution of categories/attributes of genes within each group. The categories/attributes of the genes are to be determined by the user and imported as a text (for details regarding the required format, see the [File Formats](#) section). In this analysis, hyper-geometric enrichment tests are performed, and the results can be (if requested) corrected for multiple testing using the Bonferroni correction.

To perform the analysis, select *Group Analysis >> General Enrichment Analysis >> Detect Enrichment*. The following dialog box will appear:



The following table specifies the different parameters that can be set via this dialog box:

Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Include back nodes	Include genes that are part of the module' but not included in the GE data (Relevant only if the analysis is performed on modules, detected by network based algorithm)
Load categories from	Input field for the file path, holding the gene categories/attributes.
Background set	Determines the set of genes that will be used as

	background in the analysis. Options are: all genes (of the relevant organism), original input data, filtered data or background set from file (see the Files Format section for details regarding the format of an external background set).
p-value threshold	A category/attribute will be considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than this threshold.
Multiple tests correction	Can be set to Bonferroni or None (when set to Bonferroni the corrected p-values are the ones that are compared to the threshold p-value).
Minimal overlap between category and set	The minimal number of genes from a cluster/bi-cluster expected to be categorized/attributed by an attribute in order for its enrichment to be accepted.

After the analysis is performed an enrichment analysis solution visualization tab is added to the main window. It contains general information regarding the analysis, a sort-able table holding all detected enrichments (set ID, enrichment category, p-value, etc.) and a multi-histogram panel along with a color index (mapping each color to a corresponding category). The multi-histogram panel contains one histogram for each probe/gene set/group in which enrichment has been detected. Each histogram contains a column for each significant (more frequent than would be expected by random) category. The definition of significant depends on the user's selection of threshold p-value i.e., a category is considered significantly enriched in a cluster/bicluster if its corrected p-value is lower than the preset threshold p-value.

The height of the column is proportional to the significance of this enrichment (i.e. height = $-\log(\text{raw p-value})$), and the frequency ratio (frequency in set divided by frequency in background) is written on top of the column. Upon clicking on a column, a dialog box is displayed containing the class name, corrected p-value, and a list of the genes in the cluster/bi-cluster that belong to the category. Upon clicking on one of the gene IDs in the table, a relevant web page with information regarding this gene is displayed. The display tool tip shows the cluster number, size and homogeneity.

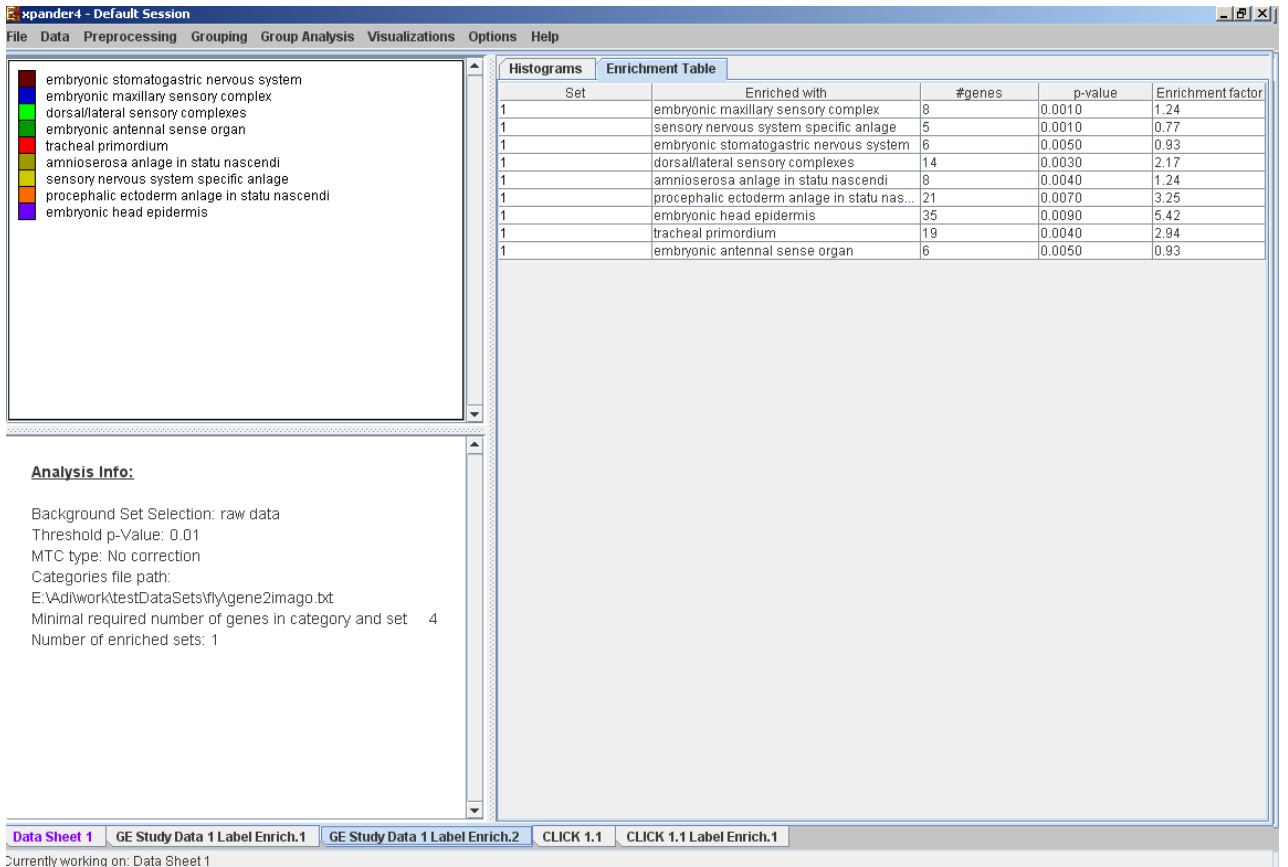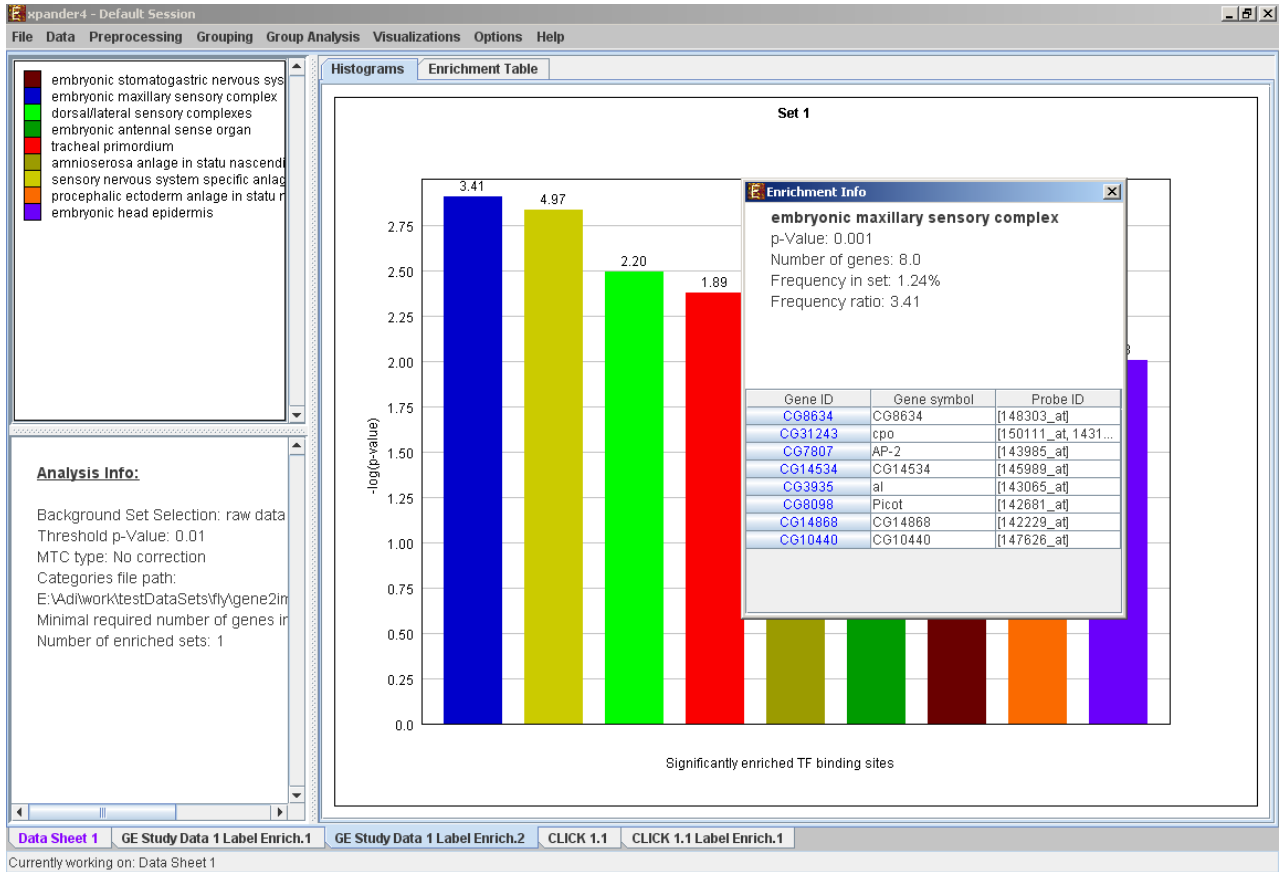


The results of this analysis can be exported to a text file (using the display menu item *Group Analysis >> General Enrichment Analysis >> Export Results*).

Network Based Group Analysis

This tool allows browsing through signaling data to view the sub-graphs that are induced by the analyzed gene groups. It also enables the user to search for statistical enrichment of these groups in highly curated signaling maps. To perform this task, Expander interfaces with the SPIKE software and database. For further information regarding the SPIKE software see the [References](#) section.

To perform the analysis on one/more of the gene groups defined in Expander (i.e. clusters, bi-clusters, modules, loaded gene sets or filtered data), select *Group Analysis >> Network >> SPIKE>>Gene Groups*. The following dialog box will appear:

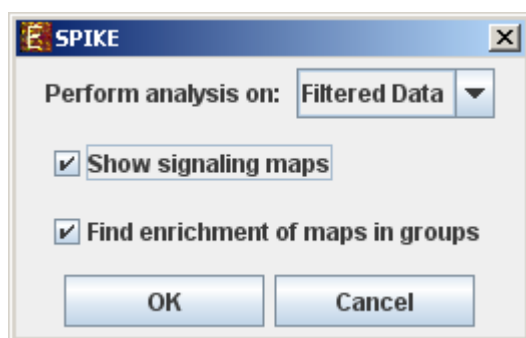

The following table specifies the different parameters that can be set via this dialog box:

Field	Description
Perform analysis on	The grouping solution on which the analysis will be performed.
Show signaling Maps	For each group display regulatory data induced by the genes included in the group.
Find enrichment of maps in groups	For each group, search for signaling maps that are enriched with genes included in the group.

Pressing OK in the dialog box will launch the SPIKE application. When operated for the first time, the launch takes a few minutes, since it has to build a local database. From this point on, please refer to page 12 in the [SPIKE user manual](#).

SPIKE can also be operated on a sub-group of genes that is derived from an existing enrichment solution in Expander. I.e. a group of genes that has a


common annotation that was found to be enriched by one of the group analysis operations. In order to operate SPIKE on such a *group*, select: *Group Analysis*>>*Network*>> *SPIKE*>>*Enrichment Derived Sets*.

Matrix Visualizations (Heat maps)

An expression matrix visualization is integrated in many of EXPANDER's displays. This visualization is similar to the red-green matrix representation of Eisen et al (1998). All it does is to render the gene-expression data on the screen in color, where green indicates under expression, and red indicates over expression. Color rendering can be configured by the user in one of the following manners: (a) by setting the range (top and bottom values) of rendered values (default values are set according to the data scale, e.g. 40-1000 for non-standardized absolute intensities data) or (b) by setting the percent of values, which are to be disregarded as extreme values from each edge (by default set to 5%). The manner of color scale configuration (i.e. (a) vs. (b)) can be set via the 'Data Matrix View' tab in the 'Display Settings' dialog box, available from *Options* >> *Settings*. The red/green coloring scheme can be changed to blue/yellow (using *Options* >> *Settings* >> *Display* >> 'Data Matrix View' tab).

A color scale appears next to the matrix (upper right side). The displayed tool tip shows the probe ID and condition title corresponding to the row and column on which the cursor is placed, and the expression value in that position. The matrix toolbar contains zoom in, zoom out, reset scale (to reset zoom factor), shorten condition title and Elongate condition title tools.

Upon selecting *Visualizations* >> *Clustered Expression Matrix*, a clustered expression matrix visualization tab is added to the main window. The probes are ordered in their original order. If a clustering solution has been previously created, its' name appears next to a radio button in the top right panel. Upon pressing this button, the order of the probes in the display changes and probe IDs are colored according to the clusters. The color index at the bottom right panel, maps each color to the index of the corresponding cluster.

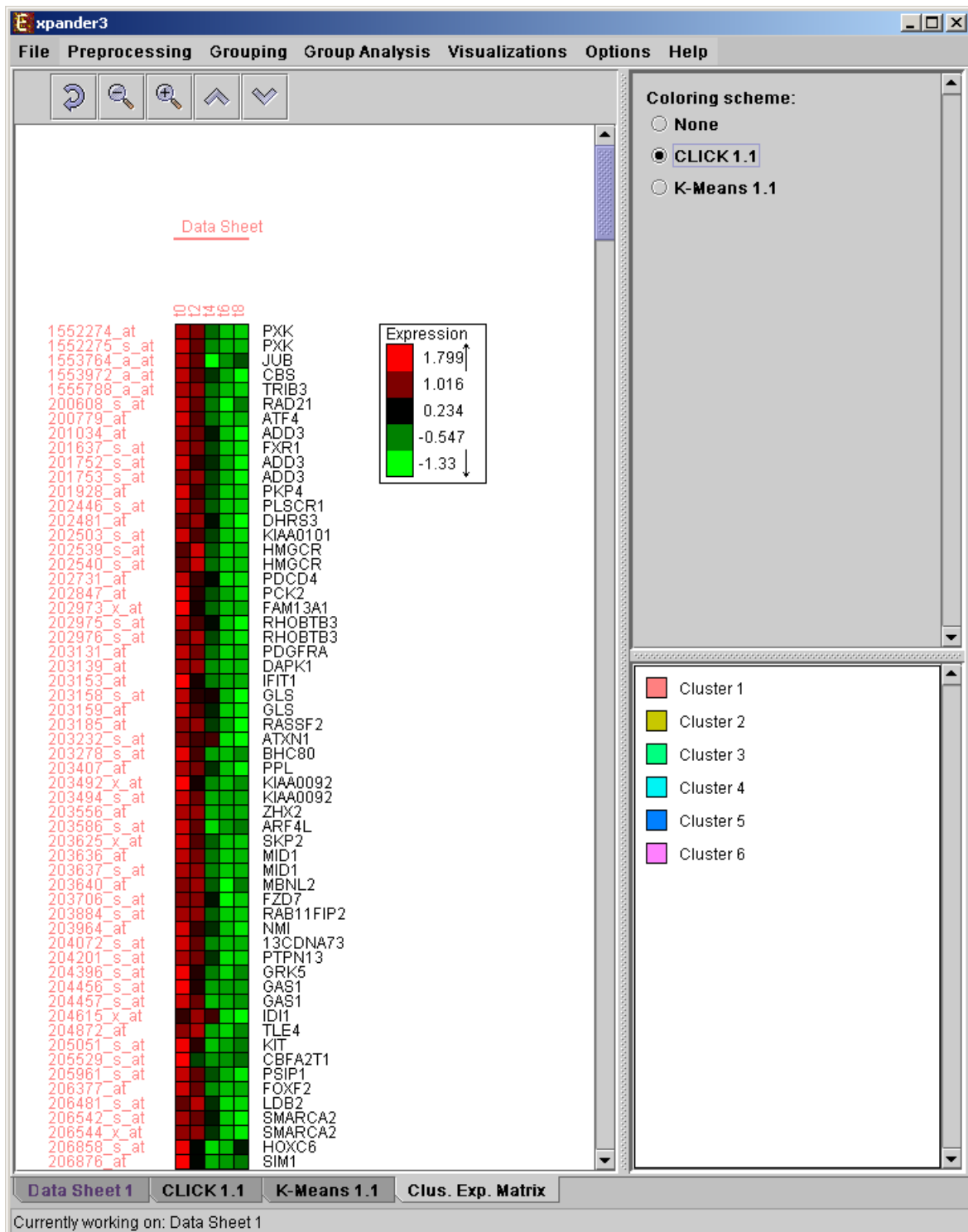


Hierarchical Clustering and Visualization

This tool uses the agglomerative algorithm to calculate a dendrogram tree for all expression patterns (probe patterns) and/or profiles (condition profiles). The type of linkage (manner in which the distance between a new node and the rest of the nodes is calculated) used in the algorithm can be set via an input dialog (for details regarding the algorithms refer to the [References](#) section). Note that it does not

generate a partition of the probes to clusters. The distance measurement used in the algorithm is $(1 - \text{Pearson Correlation})/2$.

To perform hierarchical clustering, select *Grouping >> Hierarchical Clustering*. Upon selecting this option, a dialog box appears in which the 'linkage type' parameter, used in the algorithm can be set. After pressing 'OK', the algorithm will be operated both on the probe patterns and on the condition profiles.


The resulting trees are displayed next to an expression matrix so that the probe tree appears vertically on the left and the condition tree appears horizontally above the matrix. The scale next to each tree indicates the range of distance values between vectors corresponding to the leaves. The tool tip indicates the distance value corresponding to the cursor location on the tree.

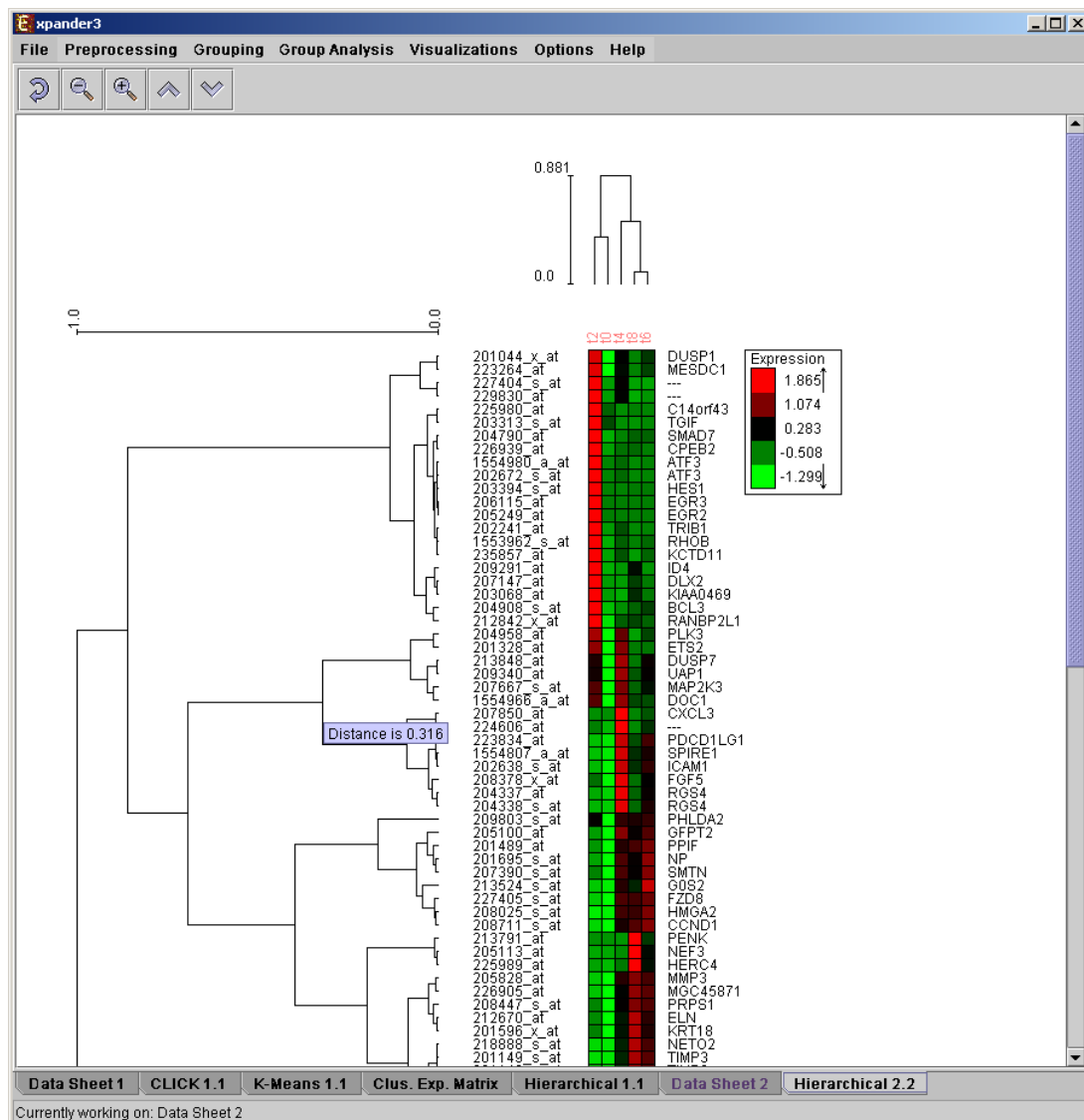

PCA Transformation

This tool transforms the original data from a k (original pattern length) to a 2 dimensional space, so that each expression vector is represented by a dot on an XY scatter chart. The transformation is based on the PCA (Principal Component Analysis) algorithm. To operate the tool, select *Visualizations >> PCA*.


If a clustering solution has been previously created, its' name appears next to a radio button in the top right panel. Upon pressing this button, the color of each dot in the display changes according to the cluster assignment of the corresponding probe. The color index at the bottom right panel, maps each color to the index of the corresponding cluster.

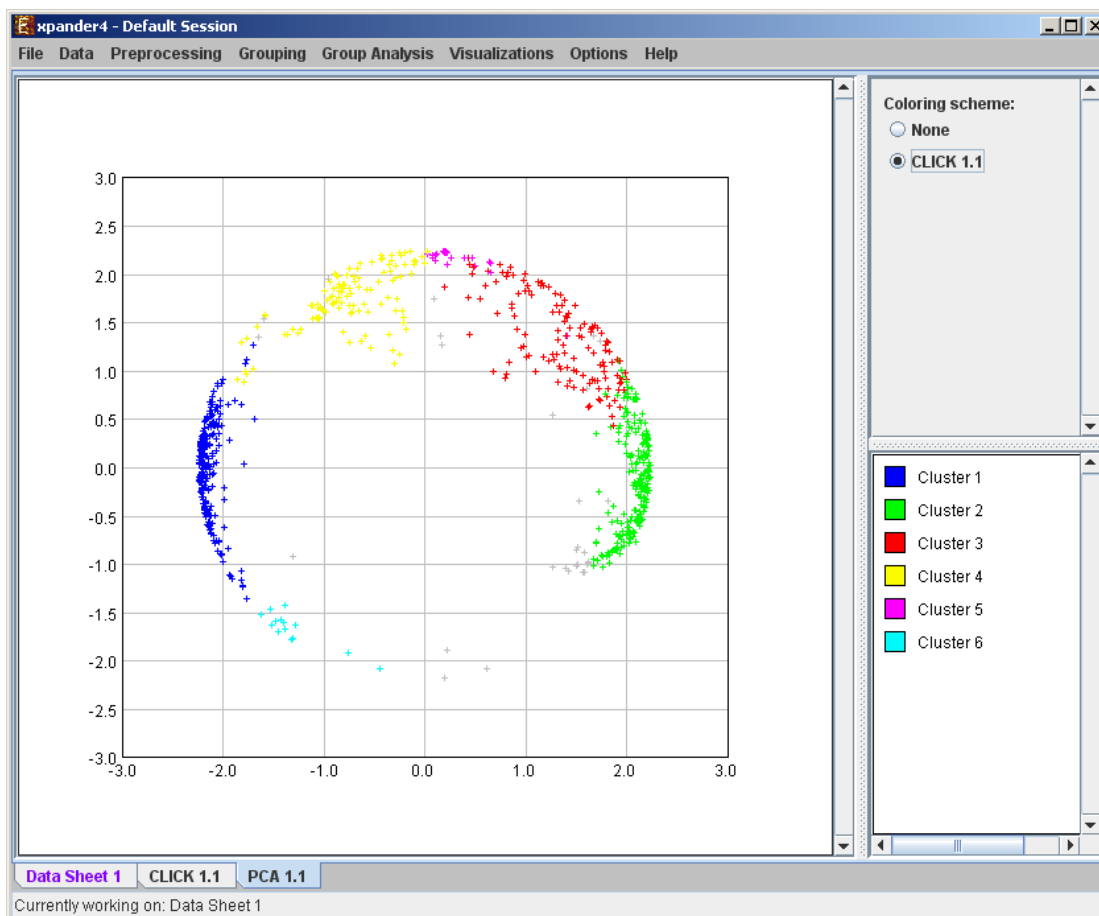

Additional Options

Searching for a gene/probe in the display


A gene can be detected in a display by selecting *Options >> Search Gene*. The following dialog box will appear:

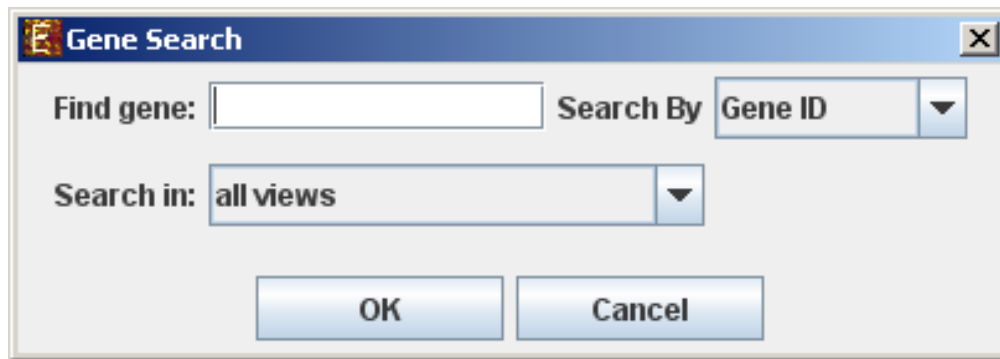

Please type the ID of the gene (can also be symbol/probe ID depending on the selection in the “Search By” combo-box) in the corresponding text box. Note that you must type the entire name or ID, not part of it. After pressing the “OK” button, a window will appear containing text describing the number of items detected in each of the searched views (number of “hits” in each view). In addition, the corresponding elements will be highlighted in all searched displays.

Defining condition subsets

You may group several conditions under a common subset name, by selecting Data >> Define Condition Subsets. This partition is used for visualization purposes. In the dialog box, select the relevant conditions, type a group name and click on the arrows.

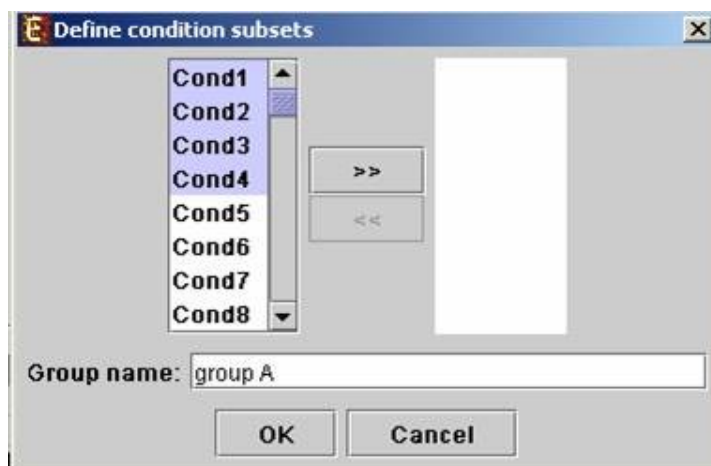

Saving and loading sessions

A set of analysis operations performed on one data set can be saved by selecting *File >> Save Session*. It can later be reloaded by selecting *File >> Load Session*. Loading a previously saved session will bring up all analysis output and visualizations that had been generated in that session, and the user will be able to continue working where he had previously stopped.

Closing views

The user can close all open views by selecting *File >> Close All*.

Closing a single view can be performed either by selecting *File>>Close* when the relevant view is selected OR by right clicking on the tab title of the relevant view and selecting *Close* from the popup menu.

Docking a view into a separate frame

Can be performed either by selecting *Options >> Dock into external frame* when the relevant view is selected OR by right clicking on the tab title of the relevant view and selecting *Dock into external frame* from the popup menu.

Upon creating the separate frame, the view will be removed from the main window. Upon closing the separate frame generated in this manner, the view will be retrieved into the main window.

Accessing the EXPANDER download page


The Expander download page can be accessed directly by selecting *Help >> Open Download Page*, while the machine is connected to the Internet.

Printing the display

Each display can be printed by selecting *File >> Print* while its tab is selected.


Exporting display into image files

Each display can be exported into image files of type .jpg, .png or .eps (post-script). This can be done by selecting *File >> Save As Image*. Upon selecting this option, a dialog box, similar to the following is displayed. In the dialog box the saved images (sections of the view), image files format, and destination directory name are input.

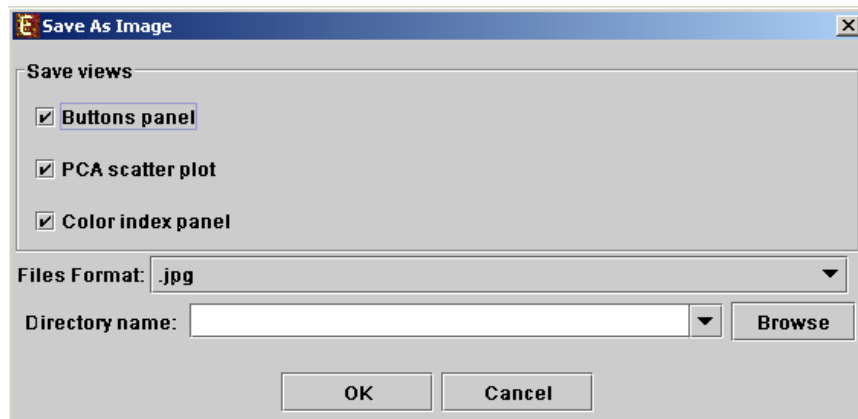

Exporting detection calls information

The detection calls info of the raw and preprocessed data can be exported into text files, by selecting *Data >> Export Detection Calls*. Upon selecting this option, the following dialog box is displayed. You may export the detection calls and also the statistics of detection calls (percent of P, M and A calls per condition), for raw data and for preprocessed data.

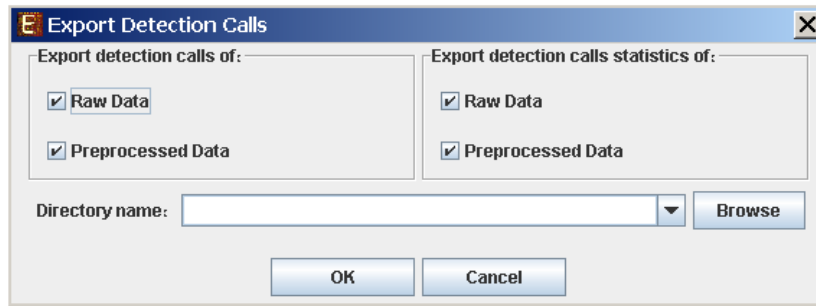

File Formats

Expression data file format:

- 1) Suffix: no limitations.
- 2) Separating token: tab delimiter.
- 3) Format:

1st line: contains a string like 'probeld' and a tab delimiter, followed by a string like 'geneSymbol' and a tab delimiter, followed by the names of all conditions separated by tab delimiters. The symbol column is optional – if the file does not contain a symbol column, please specify it in the Advanced Input Dialog box (see [Input Data](#) section).

Next lines: Each subsequent line consists of the probe ID (an identifier string that is unique to each probe in the chip), followed by a string, which represents the gene full name (if missing can be left empty by adding an additional tab delimiter), followed by its expression values (all tab delimited). If the expression file contains missing values, Expander either replaces them with a preset value (0 by default), or estimates them using the KNN (K-Nearest Neighbors) method, depending on the user selection in the data load dialog box.

*For example see files 'expressionData1.txt' and 'expressionData2.txt' in the Expander/sample_input_files/ directory.

If the data is not in the above format, it may be possible to load it using the 'Advanced' dialog box, which appears upon pressing the 'Advanced' button in the Expression Data load dialog box (see Advanced Input Dialog box in [Input Data](#) section).

Expression data with detection calls file format:

- 1) Suffix: no limitations.
- 2) Separating token: tab delimiter.
- 3) Format:

1st line: contains a string like 'probeld' and a tab delimiter, followed by a string like 'geneSymbol' and a tab delimiter, followed by the names of all conditions and detection signals columns alternately, separated by tab delimiters. Each title of condition is followed by a title of its detection column.

The symbol column is optional – if the file does not contain a symbol column, please specify it in the Advanced Input Dialog box (see [Input Data](#) section).

Next lines: Each subsequent line consists of the probe ID (an identifier string that is unique to each probe in the chip), followed by a string, which represents the gene full

name (if missing can be left empty by adding an additional tab delimiter), followed by its expression values and detection calls values, alternately (all tab delimited). Each expression value is followed by its detection value (P, M or A). If the expression file contains missing values, Expander either replaces them with a preset value (0 by default), or estimates them using the KNN (K-Nearest Neighbors) method, depending on the user selection in the data load dialog box.

*For example see files 'expressionWithDetection.txt' in the Expander/sample_input_files/ directory.

If the data is not in the above format, it may be possible to load it using the 'Advanced' dialog box, which appears upon pressing the 'Advanced' button in the Expression Data load dialog box (see Advanced Input Dialog box in [Input Data](#) section).

Gene Sets file format:

1) Suffix: no limitations

2) Format: Each line contains a gene ID, a gene symbol (optional) and the name/number of its set (separated by tabs/spaces). The gene IDs are expected to be of the same convention used in the GO annotation and TF fingerprint files. For details regarding the Gene ID convention that is used for each organism, refer to the [Supplied files](#) section.

*For example see file 'geneSetsData1.txt' under the Expander/sample_input_files/ directory (see [Sample input files](#) for more details).

Probes Filter file format:

Each line contains a single probe ID.

ID conversion file format:

1) Suffix: Currently, there are no limitations regarding the file name suffix.

2) Format: Each line contains the probe id as it appears in the data file, a tab separator and the corresponding gene ID (e.g. Entrez/Locus-Link ids for mouse and human genes and ORF codes for yeast). The second field can be left blank, indicating no conversion for that probe ID.

* It is possible that several probe IDs in the data file will be mapped to the same gene ID (e.g.: several ESTs from the same gene).

Clustering files format:

1) Suffix: no limitations.

2) Format: Each line contains the probeID, a tab separator and name/number of its cluster. The number 0 is reserved for probes that are left unclustered. The file does not have to contain all probes in the data. If a probe does not appear in the file, it is automatically set as unclustered.

*For example see file 'expressionData1Clustering.sol' (a clustering solution for the data file 'expressionData1.txt') under the Expander/sample_input_files/ directory (see [Sample Input Files](#) section for more details).

Biclustering files format:

- 1) Suffix: ``.bic``.
- 2) Format: the file is composed of two parts, presented here.

Part 1 presents a summary of the biclusters found.

- It begins with the string: ``[Bick]`` in the first line.
- Following lines contain the bicluster's id followed by its' score, separated by a tab delimiter (a line for each bicluster).

Part 2 presents the probesets and the conditions contained in each bicluster.

- It begins with the string: ``[Bicd]`` in the first line.
- Following lines contain the bicluster id, type of element ('0' for condition, '1' for probe) and element id (name of condition or probe ID), separated by tab delimiters.

Background set files format:

- 1) Suffix: no limitation.
- 2) Format: each line should contain one gene ID. The gene IDs are expected to be of the same convention used in the annotation and TF fingerprint files for the organism you are working on (please refer to the [Supplied Files](#) section).

Gene annotations/categories files format (for the general enrichment analysis):

- 1) Suffix: no limitation.
- 2) Format: each line should contain one gene ID and an annotation/category name separated by a tab delimiter. The gene IDs are expected to be of the same convention used in the annotation and TF fingerprint files for the organism you are working on (please refer to the [Supplied Files](#) section).

Sample Input Files

Several sample files are provided under `Expander/sample_input_files/`. These files include:

expressionData1.txt – A gene expression data file that was generated using the cDNA microarray technology. This is a partial dataset extracted from a yeast cell cycle dataset generated by Spellman et al 1998 (see the [References](#) section). Gene identifiers in this set are yeast ORFs, which are the same identifiers used in the annotation and TF fingerprint files that are supplied with Expander. Thus, no conversion file is required.

ExpressionData2.txt – A gene expression data file that was generated in the Affymetrix technology. This dataset was generated in an experiment that was

conducted in our laboratory on human cells, and has not yet been published. Affymetrix chips of type HG-Focus were used for this experiment and thus, the HG-Focus conversion file is required for the analysis (can be downloaded from the download page).

ExpressionData3.txt – taken from Murray JI, Whitfield ML, Trinklein ND, Myers RM, Brown PO, Botstein D: Diverse and specific gene expression responses to stresses in cultured human cells. Mol Biol Cell 2004, 15:2361-2374. A corresponding conversion file (from clone-Ids to LL-Ids) is available at the same directory under the name Data3Conversion.txt.

expressionWithDetection.txt – A gene expression data file with detection calls, that was generated in the Affymetrix technology. This dataset was generated in an experiment that was conducted in our laboratory on human cells. Affymetrix GeneChip HGU133 Plus 2.0 arrays were used for this experiment.

expressionData1Clustering.sol – A clustering solution that was generated by Expander for the dataset in 'expressionData1.txt'.

geneSetsData1.txt – Contains sets of human genes (in Entrez/Locus-Link Ids).

Data3Conversion.txt - A conversion file for expressionData3.txt.

Supplied Files

Annotation and TF-fingerprint files

Annotation and TF-fingerprint files are currently supplied with Expander for yeast, human, mouse, rat, fly, zebrafish, c-elegans and chicken. Annotation files are supplied also for E. coli and Listeria. These files are updated on a regular basis.

The following conventional gene Ids are used in the annotation and TF-fingerprint files:

Organism	ID type
Human	Entrez / Locus Link
Mouse	Entrez / Locus Link
Rat	Entrez / Locus Link
Yeast	ORF codes
Fly	FlyBase
c-elegans	WormBase
Arabidopsis	AGI Ids
Zebrafish	Ensembl
Chicken	Ensembl
E. coli*	Uniprot

Listeria*	Entrez / Locus Link
S. pombe	Entrez / Locus Link

* Only GO annotation files (no TF fingerprints data)

Gene ontology and annotation files

The gene ontology and annotation files supplied with Expander are based on data that was downloaded from the [GO website](#).

The following table specifies the database from which these files are derived for each of the organisms that we support, and the date of download from GO site:

Organism	Extracted from	Updated on
Human	GOA@EBI	May 2007
Mouse	MGI (Mouse Genome Informatics)	May 2007
Rat	GOA@EBI	October 2007
Yeast	SGD (Saccharomyces Genome Database)	October 2007
Fly	FlyBase	January 2005
c-elegans	WormBase	October 2007
Arabidopsis	TAIR (The Arabidopsis information resource)	October 2007
Zebrafish	ZFIN	May 2007
Chicken	GOA@EBI	January 2008
e.coli	EBI	April 2008
Listeria	Blast2Go	February 2009
S. pombe	NCBI	November 2008

TF fingerprint files and sequence files

The following TF fingerprint files (and sequence files) can be downloaded from our web download page (as part of the zip supplied for each organism).

Organism	Promoter sequences range	Promoter sequences origin	TF Models origin
Human	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 42	TRANSFAC (version 8.2)
Mouse	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 42	TRANSFAC (version 8.2)
Rat	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 47	TRANSFAC (version 8.2)
Yeast	600 bp upstream CDS to CDS	SGD January 2007	TRANSFAC (version 8.2) & Harbison et al. 2004 (see References)
Fly	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 43	TRANSFAC (version 8.2)
C-elegans	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 44	TRANSFAC (version 8.2)
Arabidopsis	1000 bp upstream the TSS to 200 bp downstream the TSS.	TAIR: The Arabidopsis Information Resource (December '06)	TRANSFAC (version 8.2)
Zebrafish	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 44	TRANSFAC (version 8.2)
Chicken	3000 bp upstream the TSS to 200 bp downstream the TSS.	Ensembl release 43	TRANSFAC (version 8.2)
S. pombe	600 bp upstream CDS to CDS	Sanger gene DB October 2008	TRANSFAC (version 8.2) & Harbison et al. 2004

miRNA target scan files:

The following miRNA target scan files can be downloaded from our web download page (as part of the zip supplied for each organism).

Organism	Extracted from	Updated on
Human	TargetScan website	TargetScan5
Mouse	TargetScan website	TargetScan5
Fly	TargetScan website	TargetScan5
c-elegans	TargetScan website	TargetScan5

Gene positions data files:

The following gene position files can be downloaded from our web download page (as part of the zip supplied for each organism). It is used to display the gene chromosomal positions.

Organism	Extracted from	Updated on
Human	UCSC genome browser	January 2009
Mouse	UCSC genome browser	January 2009
Rat	UCSC genome browser	January 2009
Fly	UCSC genome browser	January 2009
c-elegans	UCSC genome browser	January 2009
Zebrafish	UCSC genome browser	January 2009

Gene ID conversion files:

Gene ID conversion files for many of the Affymetrix chips can be downloaded from the Expander download page. The files map each Affymetrix Id into the corresponding gene Id. Conversion files are generated and added to the download page according to user requests. If you can't find the file you need here, please look it up in the download page, and [contact](#) us if it's not there.

Organism	Chip name
Human	HG-Focus
Human	HGU1332
Human	HG-U95E
Human	HG-U133A
Human	HT_HG-U133A
Human	HG-U133Plus2
Human	Hu-35KsubB
Human	HuGene-1_0-ST
Mouse	MGU74Av2
Mouse	MGU430_2
Mouse	MG430A2
Mouse	MoGene-1_0-ST
Rat	RGU34A
Rat	Rat230_2
Rat	Agilent
C-elegans	C. elegans Genome Chip
Arabidopsis	ATH1
Zebra-Fish	GeneChip Zebrafish Genome Array
Chicken	Affymetrix Chicken Genome Chip

E. coli	Affymetrix E. coli Antisense Genome Array
E. coli	Affymetrix E. coli Genome 2.0 Array

Network files :

Organism	File name	Network origin
Human	Expander.hsa.RualNature05.sif	Towards a proteome-scale map of the human protein-protein interaction network by Rual JF et al. <i>Nature</i> . 437(7062):1173-8 (2005)
Human	Expander.hsa.IntAct.sif	IntAct database (http://www.ebi.ac.uk/intact/)
Mouse	Expander.mmu.IntAct.sif	IntAct database (http://www.ebi.ac.uk/intact/)
Rat	Expander.rno.IntAct.sif	IntAct database (http://www.ebi.ac.uk/intact/)
Worm	Expander.cel.SimonisNatMethods08.sif	Empirically controlled mapping of the <i>Caenorhabditis elegans</i> protein-protein interactome network by Simonis N. et al. <i>Nature Methods</i> 6, 47 - 54 (2009)
Fly	Expander.dme.Droid.sif	Droid database (http://www.droidb.org/)
Yeast	Expander.sce.United.sif	<ol style="list-style-type: none"> 1. High-Quality Binary Protein Interaction Map of the Yeast Interactome Network by Yu et al. <i>Science</i> 322(5898):104 – 110 (2008) 2. Comprehensive curation and analysis of global interaction networks in <i>Saccharomyces cerevisiae</i> by Reguly et al. <i>Journal of Biology</i> 5(4):11 (2006) 3. Toward a comprehensive atlas of the physical interactome of <i>Saccharomyces cerevisiae</i> by Collins SR et al. <i>Molecular Cell Proteomics</i> 6(3):439-50 (2007)
Arabidopsis	Expander.ath.TAIR.sif	TAIR database (http://www.arabidopsis.org/)
E. coli	Expander.eco.Arifuzzaman06.txt	Large-scale identification of protein-protein interaction of <i>Escherichia coli</i> K-12 by Arifuzzaman M et al. <i>Genome Research</i> 16(5):686-91. (2006)

Settings

The Settings are accessible from the *Options* menu, and contain *Display* settings *External applications* settings.

The **Display** dialog box contains the following tabs:

Clustering Results View – Contains check boxes that configure the following parameters:

- A common Y-axis scale for all cluster patterns (vs. cluster specific)
- Visible x axis
- Connect all points in a pattern
- Display similarity matrix for probes (using Pearson correlation)
- Display similarity matrix for conditions (using Pearson correlation)

Group Analysis Results View – Contains a check box that configures whether the Y-axis scale of all histograms is common OR cluster specific.

Data Matrix View – Allows selection between:

- Range control and extreme values control when rendering expression matrix values.
- The red/green coloring scheme can be changed to blue/yellow.

The **External applications** dialog box allows specification of the location of the R executable (required for CEL files loading). In Windows, R.exe file is likely to be located in the 'bin' folder of R software. In Linux, you may type 'which R' in the command line to find R path.

R External Application

The CEL file preprocessing and the newly added SAM filter utilities require the pre-installation of one of the recent versions of R, a free software environment for statistical computing and graphics. R can be installed from: <http://cran.r-project.org/>.

After installing R, please do the following to install the Bioconductor “affy” package, “aroma.affymetrix” package, “gcrma” package and the “samr” package:

1. Run R.
2. In the R frame\window type the text:
source("http://bioconductor.org/biocLite.R")
3. Press ‘Enter’.
4. In the R frame\window type the text: **biocLite("affy")**
5. Press ‘Enter’.

To install the 'samr' package:

6. In the R frame\window type the text: **install.packages("samr")**

7. Press 'Enter'.

To install the 'aroma.affymetrix' package:

8. In the R frame\window type the text:
source("http://www.braju.com/R/hbLite.R");

9. Press 'Enter'.

10. In the R frame\window type the text: **hbInstall("aroma.affymetrix");**

To install the 'gcrma' package:

11. In the R frame\window type the text:
source("http://bioconductor.org/biocLite.R")

12. Press 'Enter'. In the R frame\window type the text:

13. biocLite("gcrma")

14. Press 'Enter'.

You may install only one of the packages, depending on what you wish to use (to install only "samr", follow instructions number 1, 6 and 7).

If Expander cannot find your R software, a window will pop, asking you to specify its location. Please browse to the location of your R software. In Windows, R.exe file is likely to be located in the 'bin' folder of R software. In Linux, you may type 'which R' in the command line to find R path. If you have a few versions of R installed, please make sure to point Expander to a version in which the necessary packages have been installed.

You may also specify R location from the menu: *Options >> Settings >> External applications.*

To use R utilities, please make sure there are no white spaces in the path of Expander directory (or the CEL files directory, if loading CEL files). For example, if the name of the Expander folder is 'Expander 5', change it to 'Expander_5'. If Expander is under "Program Files" it should be moved to another location, because of the space between "Program" and "Files". The R software does not cope well with spaces in the path.

Also, please make sure to have 'write' permission to the Expander\Rscripts directory. If you are loading CEL files, check also that you have 'write' permission to the *Files location* which you specified in the 'Load CEL Files' dialog box.

FAQ

Linux/Unix problems

[Click clustering, Samba bi-clustering, Tango and Prima Promoter analysis algorithms fail when running on Linux/Unix.](#)

CEL Files Loading Problems

[How do I install R and the Bioconductor "affy" package?](#)

[Loading of CEL files or performing SAM filter continue for ever](#)

[Loading of CEL files fail.](#)

Clustering

[When I try to run Biclustering on my data I get a failure notice.](#)

[How can I save the clustering expression patterns charts?](#)

Grouping Analysis (functional and promoter analysis)

[When I run Functional Analysis, Expander gets stuck.](#)

[When I load a session with that contains Functional Analysis results, Expander gets stuck,](#)

[When I try to run Promoter Analysis no values appear in the Fingerprints file field of the input dialog box.](#)

[When I try to run the promoter/functional analysis, I get a failure message box.](#)

[Promoter/Functional analysis produces no results \(the resulting view is empty\).](#)

[How can I save the bar charts produced by Expander, displaying the enrichments?](#)

[Why do certain Transcription Factors have a few accession numbers or/and a few gene IDs?](#)

Others

[Can I run Expander on Mac OS?](#)

Click clustering, Samba bi-clustering, Tango and Prima Promoter analysis algorithms fail when running on Linux/Unix.

Answer: Make sure that you have write permission in the Expander directory, and execution permissions on the files: click.exe, samba.exe, annot_sets.exe and analyzeFingerprints.exe, which are under the Expander directory. If the problem still occurs, open the file expanderLog.txt and search for the text: "libstdc++.so.5". If this text appears (along with a message indicating it has not been found), please contact your system administrator and report this problem (this is a system problem). If you do not have a system administrator, and fail to install this library, please contact us (expander@cs.tau.ac.il) and we will try to assist.

[How do I install R and the Bioconductor "affy" package?](#)

Answer: please refer to [R External Application](#) section.

[Loading of CEL files or performing SAM filter continue forever.](#)

Answer: If the operation continues forever (the 'processing, please wait' window is displayed), please check if there is a folder with a space in its name somewhere in the path of Expander (or the CEL files) directory. For example, if the name of the Expander folder is 'Expander 4', change it to 'Expander_4'. The R software used for preprocessing CEL files has a problem dealing with spaces in the path. If this is the problem, then in the expanderLog.txt file (in your Expander directory) there should be a message about arguments being ignored.

[Loading of CEL files fail.](#)

Answer: If you are working with Expander version 4, please update to 4.0.1 or higher. Make sure you have R along with the Bioconductor "affy" / "aroma.affymetrix" packages installed in the R version which is specified in the settings "External Applications" tab (from the menu select *Options >> Settings >> External applications*). If R location is not defined in the settings, please define it (In Windows, R.exe file is likely to be located in the 'bin' folder of R software. In Linux, you may type 'which R' in the command line to find R path). If loading of CEL files still fails, please make sure that the *Files location* which you specified in the 'Load CEL Files' dialog box, is a folder which contains CEL files and that you have write permission to that folder.

When I try to run Biclustering on my data I get a failure notice.

Answer: Make sure that the 'Use option files of type' field in the SAMBA input dialog box is not empty (if it is, please re-download Expander). Also make sure that the following files exist in your Expander directory: *ibic.opt*, *samba.exe*.

How can I save the clustering expression patterns charts?

When the clustering results tab is open, please go to *File >> Save As Image*.

When I run Functional Analysis, Expander gets stuck.

Answer: If you are working with Expander version 4.0 or 4.0.1, please update to a higher version (4.0.2 and on).

When I load a session with that contains Functional Analysis results, Expander gets stuck.

Answer: If you are working with Expander version 4.0 or 4.0.1 (and the session was created with a version < 4.0), please update to a higher version (4.0.2 and on).

When I try to run Promoter Analysis no values appear in the Fingerprints file field of the input dialog box.

Answer: Fingerprint files are not placed in the right directory. Fingerprint files should be placed under the 'TF_fingerprints' directory that is under the Expander/organism/<org name> directory. For example, the human FP file should be placed under: *.../Expander/organisms/human/TF_fingerprints/*. When downloading the organism specific data zip, it should be extracted into the Expander/organisms/ directory. This will automatically put them in the right place.

Group analysis (enrichment analysis) leads to a failure message box.

Answer: Errors while running group (enrichment) analysis can be caused by the following problems:

- a) Organism specific data (Fingerprint\annotation files) is not in the right directory. The organism specific data zip should be extracted into the Expander/organisms/ directory. You may download the relevant data by selecting from the menu: *Help >> Download Data for Organism*.

- b) Data contains elements that do not appear in the background set (this is only relevant when the background set is loaded from an external file).

Group analysis (enrichment analysis) produces no results (the resulting view is empty).

Answer: This can be caused by one of the following:

a) You are using the wrong conversion file or a conversion file that is not in the right format or does not map the probes to the expected type of gene IDs. The conversion file maps each probe ID in your data file to a gene ID that is used in the group analysis. A conversion file is required when the probe IDs in your data file do not match the ones in the enrichment files (for example annotation and TF_fingerprint files that we supply).

b) You did not set the organism field in the input dialog to the organism type of your data.

c) You are trying to analyze only one set (e.g. the filtered data set) which you are using also as background (in this case the analysis has no meaning since it is trying to detect enrichments in the cluster/bicluster in comparison to the background set).

d) You set the threshold p-value to be too strict (low).

e) Biological reason i.e., there is nothing to report regarding this specific clustering/biclustering solution or this gene sets data.

How can I save the bar charts produced by Expander, displaying the enrichments?

When the results tab is open, please go to *File >> Save As Image*.

Why do certain Transcription Factors have a few accession numbers or/and a few gene IDs?

Answer: The transcription factors (TFs) found enriched by Prima are presented in the following way:

Accession Num. in TRANSFAC DB [TF name]. For example, M00287[NF-Y]

It is possible that a TF will have a few accession numbers in TRANSFAC, which represent different PWMs (position weight matrices specify the probability for observing each nucleotide at each position of the binding site, based on a set of empirically validated binding sites of the respective TF).

It is also possible that a TF will have a few Entrez gene IDs, since a TF may be composed of a few proteins. For example, NF-Y is a trimer, composed of 3 subunits.

Can I run Expander on Mac OS?

Expander is not designed for Mac OS. You can probably use it partially – without running its features that require the execution of exe files (CLICK, SAMBA, TANGO and PRIMA). The exe files are only suitable for Windows and Linux / Unix.

This section will be updated as we get user feedbacks and problems.

Please refer all questions/comments to Expander@cs.tau.ac.il.

Copyrights Information

Copyrights © Tel-Aviv University, Israel (2003).

This product uses the FreeHEP Java Library, which is distributed under the [LGPL](#) license. FreeHEP copyright holders: CERN, Geneva, Switzerland SLAC, Stanford, California, U.S.A. University of California Santa Cruz, U.S.A.

This product uses the XStream Java Library, which is distributed under the BSD license (see BSD.txt). Copyright holders: (c) Joe Walnes 2003-2005

This product uses the Caryoscope java component, which is distributed under the MIT license (see MIT_Caryoscope.txt). Copyright holders: Copyright 2003-2004, Ihab A.B. Awad; Copyright 2006, Anjalee Sujanani; Stanford University.

A portion of the user interface code is due to Sun Microsystems, Inc. Copyright 1994-2004 Sun Microsystems, Inc. All Rights Reserved. The following license rules apply to that portion:

'Neither the name of Sun Microsystems, Inc. or the names of contributors may be used to endorse or promote products derived from this software without specific prior written permission.

This software is provided 'AS IS,' without a warranty of any kind. ALL EXPRESS OR IMPLIED CONDITIONS, REPRESENTATIONS AND WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR NON-INFRINGEMENT, ARE HEREBY EXCLUDED. SUN MICROSYSTEMS, INC. ('SUN') AND ITS LICENSORS SHALL NOT BE LIABLE FOR ANY DAMAGES SUFFERED BY LICENSEE AS A RESULT OF USING, MODIFYING OR DISTRIBUTING THIS SOFTWARE OR ITS DERIVATIVES. IN NO EVENT WILL SUN OR ITS LICENSORS BE LIABLE FOR ANY LOST REVENUE, PROFIT OR DATA, OR FOR DIRECT, INDIRECT, SPECIAL, CONSEQUENTIAL, INCIDENTAL OR PUNITIVE DAMAGES, HOWEVER CAUSED AND REGARDLESS OF THE THEORY OF LIABILITY, ARISING OUT OF THE USE OF OR INABILITY TO USE THIS SOFTWARE, EVEN IF SUN HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.

You acknowledge that this software is not designed, licensed or intended for use in the design, construction, operation or maintenance of any nuclear facility.'

References

R: R Development Core Team (2008). R: A language and environment for statistical computing. R Foundation for Statistical Computing, Vienna, Austria. ISBN 3-900051-07-0, URL <http://www.R-project.org>

affy package for CEL files preprocessing: Gautier, L., Cope, L., Bolstad, B. M., and Irizarry, R. A. 2004. affy---analysis of Affymetrix GeneChip data at the probe level. *Bioinformatics* 20, 3 (Feb. 2004), 307-315

aroma.affymetrix package for CEL files preprocessing:

H. Bengtsson; K. Simpson; J. Bullard; K. Hansen. aroma.affymetrix: A generic framework in R for analyzing small to very large Affymetrix data sets in bounded memory, Tech Report #745, Department of Statistics, University of California, Berkeley, February 2008

H. Bengtsson, R. Irizarry, B. Carvalho, & T.P. Speed. Estimation and assessment of raw copy numbers at the single locus level, *Bioinformatics*, 2008

gcrma R package: Jean(ZHIJIN) Wu and Rafael Irizarry with contributions from James MacDonald Jeff Gentry (). gcrma: Background Adjustment Using Sequence Information. R package version 2.14.1.

Quantile normalization: Bolstad, B. M. Irizarry, R. A. Astrand, M. and Speed, T. P. A Comparison of Normalization Methods for High Density Oligonucleotide Array Data Based on Variance and Bias. *Bioinformatics* 19(2):185-193, 2003

Non-linear baseline normalization: Schadt, E., C. Li, B. Eliss, and W. H. Wong. Feature extraction and normalization algorithms for high-density oligonucleotide gene expression array data. *J. Cell. Biochem.* 84(S37),120–125, 2002

SAM (Significance Analysis of Microarray):

V. Tusher., R. Tibshirani., and G. Chu. Significance analysis of microarrays applied to the ionizing radiation response. *PNAS*, 98: 5116-5121, 2001

R. Tibshirani, G. Chu, T. Hastie and Balasubramanian Narasimhan (). samr: SAM: Significance Analysis of Microarrays. R package version 1.26. <http://www-stat.stanford.edu/~tibs/SAM>

K-Means clustering algorithm: Tavazoie, S., Hughes, J. D., Campbell, M. J., Cho, R. J., and Church, G. M. Systematic determination of genetic network architecture. *Nat Genet*, 22: 281-285, 1999

SOM clustering algorithm: Tamayo, P., Slonim, D., Mesirov, J., Zhu, Q., Kitareewan, S., Dmitrovsky, E., Lander, E. S., and Golub, T. R. Interpreting patterns of gene expression with self-organizing maps: methods and application to hematopoietic differentiation. *Proc Natl Acad Sci U S A*, 96, 2907-2912, 1999

CLICK clustering algorithm: Sharan, R. and Shamir, R. CLICK: a clustering algorithm with applications to gene expression analysis. *Proc Int Conf Intell Syst Mol Biol* 8, 307-16, 2000

SAMBA biclustering algorithm: Tanay, A. Sharan, R. and Shamir, R. Discovering statistically significant biclusters in gene expression data. *Bioinformatics*, 18(1), 136-144, 2002

Matisse network grouping: Ulitsky, I. and Shamir, R. MATISSE: Identification of functional modules using network topology and high-throughput data. *BMC Systems Biology*, vol 1, No. 8 (2007)

Context Scores for miRNA enrichment analysis: Grimson, A., Kai-How Farh, F. K Johnston, W., Garrett-Engle, F., P Lim, L., P Bartel, D. MicroRNA Targeting

Specificity in Mammals: Determinants beyond Seed Pairing. *Molecular Cell*, 27:91-105 (2007)

PRIMA algorithm: Elkon, R., Linhart, C. Sharan, R. Samir, R. and Shiloh, Y. Genome-Wide In Silico Identification of Transcriptional Regulators Controlling the Cell Cycle in Human Cells. *Genome Research*, Vol. 13(5), pp. 773-780, 2003.

Spike software and DB: R. Elkon, R. Vesterman, N. Amit, I. Ulitsky, I. Zohar, M. Weisz, G. Mass, N. Orlev, G. Sternberg, R. Blekhman, J. Assa, Y. Shiloh and R. Shamir. SPIKE - a database, visualization and analysis tool of cellular signaling pathways. *BMC Bioinformatics* 2008, **9:11**

Spike home page: <http://www.cs.tau.ac.il/~spike/>

Agglomerative algorithm for hierarchical clustering: Eisen, M. B., Spellman, P. T. et al. Cluster analysis and display of genome-wide expression patterns. *Proc Natl Acad Sci U S A* 95(25), 14863-8, 1998

TF binding site profiles that were used to generate the supplied yeast TF fingerprint files: Harbison, C.T., D.B. Gordon, T.I. Lee, N.J. Rinaldi, K.D. Macisaac, T.W. Danford, N.M. Hannett, J.B. Tagne, D.B. Reynolds, J. Yoo, E.G. Jennings, J. Zeitlinger, D.K. Pokholok, M. Kellis, P.A. Rolfe, K.T. Takusagawa, E.S. Lander, D.K. Gifford, E. Fraenkel, and R.A. Young. Transcriptional regulatory code of a eukaryotic genome. *Nature* 431: 99-104, 2004

expressionData1.txt sample input file: Spellman, P. T., Sherlock, G., et al. Comprehensive identification of cell cycle-regulated genes of the yeast *Saccharomyces cerevisiae* by microarray hybridization. *Mol Biol Cell* 9(12), 3273-97, 1998